

vegetables and fruits by modifying the official Multi-residue method II (Fig. 1b) implemented in Japan. The sample preparation procedure was carefully optimized to develop a highly efficient method for the removal of co-extracted matrix components. The newly developed method was applied for the simultaneous determination of 73 multi-class acidic pesticides in vegetables and fruits.

## Experimental

### Reagents and materials

**Solvents and chemicals.** Pesticide-analysis-grade acetone, acetonitrile, hexane, methanol, and toluene were obtained from Kanto Chemical (Tokyo, Japan). Pesticide-analysis-grade sodium chloride and analytical-grade ammonium sulfate, formic acid, and 0.1 mol L<sup>-1</sup> hydrochloric acid were purchased from Wako Pure Chemical (Osaka, Japan). LC-MS-grade methanol and water were obtained from Kanto Chemical and used for LC-MS/MS analyses. Diatomaceous earth (Celite<sup>®</sup>, No. 545) was obtained from Wako Pure Chemical. Water used to prepare test solutions was purified in an NZJ-2DSYW distillation apparatus (Fujiwara Scientific, Tokyo, Japan).

**Analytical standards.** Pesticide standards were purchased from Dr. Ehrenstorfer (Augsburg, Germany), Hayashi Pure Chemical (Osaka, Japan), Wako Pure Chemical, and Kanto Chemical. Individual stock standard solutions (1,000 mg L<sup>-1</sup>) were prepared in acetonitrile or methanol, depending on the solubility of pesticide in each solvent. Working solutions were prepared by mixing the stock standard solutions and diluting with methanol. Calibration standard solutions were freshly prepared by diluting working standard solutions with methanol.

**Cartridge columns.** Octadecylsilyl silica gel (ODS) columns (Mega Bond Elut C18, 1,000 mg) were obtained from Agilent Technologies (Palo Alto, CA, USA), and graphitized carbon (InertSep GC, 500 mg) and silica gel (InertSep SI, 1 g) columns were obtained from GL Sciences (Tokyo, Japan).

### Apparatus

**LC-MS/MS.** An Acquity UPLC system (Waters, Milford, MA, USA) coupled with a Xevo TQ-S mass spectrometer (Waters) was used, with the following operating conditions: Inertsil ODS-4 column (length: 150 mm, i.d.: 2.1 mm, particle size: 3  $\mu$ m; GL Sciences); mobile phases, 5 mmol L<sup>-1</sup> ammonium acetate in water (Solvent A) and 5 mmol L<sup>-1</sup> ammonium acetate in methanol (Solvent B); solvent gradient, 10% Solvent B at 0 min, 95% Solvent B at 20 min, 95% Solvent B at 30 min, 100% Solvent B at

30.1 min, 100% Solvent B at 40 min, 10% Solvent B at 40.1 min; flow rate, 0.2 mL min<sup>-1</sup>; column temperature, 40°C; injection volume, 5  $\mu$ L; ionization mode, electrospray ionization (ESI); capillary voltage, 3 kV (ESI (+)), 1 kV (ESI (-)); source temperature, 150°C; desolvation temperature, 500°C; desolvation gas, nitrogen at 1,000 L h<sup>-1</sup>; cone gas, nitrogen at 150 L h<sup>-1</sup>; and collision gas, argon at 0.15 mL min<sup>-1</sup>. The selective reaction monitoring (SRM) transitions and retention times are summarized in Table A1.

**Laboratory apparatus.** A homogenizer (Polytron PT 10-35 GT; Kinematica, Lucerne, Switzerland), food processor (Grindomix GM 200; Retsch, Haan, Germany), rotary evaporator (N-1000/NVC-2100; Tokyo Rikakikai, Tokyo, Japan), electric shaker (SR-2W; Taitec, Saitama, Japan), and vacuum pump (APN-215MV-1; Iwaki, Tokyo, Japan) were connected to a Kiriya funnel containing filter paper No. 5B (Kiriya Glass Works, Tokyo, Japan).

### Sample preparation

**Foods containing low amounts of pigments.** A 20.0-g sample was weighed in a 250-mL glass tube, and was extracted with 40 mL of acetonitrile and 10 mL of 0.1 mol L<sup>-1</sup> hydrochloric acid by using a homogenizer. The homogenate was filtered under vacuum through Celite<sup>®</sup> filter aid, and the residue was re-homogenized with 20 mL of acetonitrile and filtered again. The volume of the combined extract was made up to 100 mL with acetonitrile.

A 10-mL aliquot of the extract was added to a 50-mL PTFE centrifuge tube containing 5 g of sodium chloride and shaken vigorously for 10 min. The extract was then centrifuged at 3,000 rpm for 5 min. The resultant acetonitrile layer was loaded on an ODS column preconditioned with 10 mL of acetonitrile and eluted with an additional 5 mL of acetonitrile. The combined eluate was concentrated to approximately 1 mL with a rotary evaporator below 40°C and evaporated to dryness under a stream of nitrogen.

The residue was subsequently re-dissolved in 2 mL of hexane/acetone/formic acid (100:100:1 v/v/v) and loaded on a silica gel column (1 g), eluting with an additional 18 mL of hexane/acetone/formic acid (100:100:1 v/v/v). The combined eluate was concentrated to approximately 1 mL with a rotary evaporator below 40°C, evaporated to dryness under a stream of nitrogen, and the residue was re-dissolved in 4 mL of methanol.

**For foods containing high amounts of pigments (spinach).** Samples were prepared using the same procedure as mentioned above. The residue obtained after silica gel column cleanup was re-dissolved in 2 mL of acetonitrile/toluene/water (30:10:1 v/v/v), and subsequently loaded onto a

graphitized carbon column (500 mg) and eluted with an additional 23 mL of acetonitrile/toluene/water (30:10:1 v/v/v). The combined eluate was concentrated to approximately 1 mL with a rotary evaporator below 40°C and evaporated to dryness under a stream of nitrogen; the obtained residue was re-dissolved in 4 mL of methanol.

#### Recovery tests

The recovery tests were performed for cabbage, spinach, potato, eggplant, orange, and apple, with five replicates for all matrices at a spiking level of 0.01 mg kg<sup>-1</sup>. A 1-mL aliquot of working standard solution was added to the samples and were allowed to stand for 30 min before extraction. Quantifications were carried out using six points (12.5, 25, 50, 75, 100, 150%) from external calibration using solvent standards.

*Preparation of matrix-matched standards.* Matrix-matched standards were prepared as follows: Blank extract solutions of 100 µL were evaporated to dryness under a stream of nitrogen, and the residue was re-dissolved in 100 µL of 0.005 µg mL<sup>-1</sup> standard mixture (0.0005 µg mL<sup>-1</sup> for warfarin) in methanol.

## Results and discussion

#### Selection of pesticides and optimization of LC-MS/MS conditions

In this study, a total of 73 LC-amenable acidic pesticides, mostly herbicides with a wide pK<sub>a</sub> range,<sup>[2-6]</sup> were selected (Table A1). These pesticides are retained on PSA or NH<sub>2</sub> sorbent when using acetonitrile/toluene (3:1 v/v) as elution solvent, and could not be determined by the official Japanese Multi-residue method I (Fig. 1a).

Optimizations of LC-MS/MS conditions were performed by flow injection analysis of individual standards. The precursor and product ions were optimized by varying cone voltage in the range of 10–80 V and collision energy in the range of 5–45 eV in both ESI(+) and ESI(–) modes. The most intense transition was used for quantification, and the second transition was used for confirmation (Table A1).

Figure 2 shows representative chromatograms of standards, blank apple sample, and fortified sample. Satisfactory separations and peak shapes were achieved for most of the tested pesticides by gradient elution using methanol and water containing 5 mmol L<sup>-1</sup> of ammonium acetate, except for highly polar pesticides with short retention time (i.e. mesotrione, picloram, and trinexapac), which gave tailing peaks.

#### Optimization of sample preparations

To develop a reliable and rugged Multi-residue method that can sufficiently extract acidic pesticides having a wide range of pK<sub>a</sub> values from foods, regardless of their pH, the pH of extraction should be controlled.<sup>[12]</sup> In addition, homogenization with a probe blender is known to accelerate extraction more significantly than shaking.<sup>[21]</sup> Thus, in this study, 0.1 mol L<sup>-1</sup> hydrochloric acid was added to acetonitrile, and extraction was performed by homogenizing with a probe blender.

To remove water from the crude extract by salting out, 5 g of sodium chloride was added to 10 mL of crude extract and shaken for 10 min. Since a mixture of hydrochloric acid and acetonitrile was used for extraction in this study, the pH values of aqueous (lower) layers after the salting-out step are as follows: cabbage, 3.3; spinach, 2.2; potato, 3.4; eggplant, 2.3; orange, 2.2; and apple, 2.1. All pesticides were partitioned into the acetonitrile layer in all the tested foods, although the recoveries of seven pesticides out of 73 were low under neutral conditions (pH 7). It seems that in this pH range most acidic pesticides exist in their neutral forms; thus, it was unnecessary to add additional acid at the salting-out step.

The presence of low-polarity matrix components during determination would be troublesome, since they affect the LC-MS/MS performance and cause matrix effects; therefore, these components should be sufficiently removed prior to LC-MS/MS analysis. Thus, in this study, ODS column cleanup was added after the salting-out step. The acetonitrile (upper) layer obtained after salting out was directly loaded on an ODS column and eluted with an additional 5 mL of acetonitrile. Low-polar matrix components, such as green pigments (chlorophylls), were retained on the column, while all the tested pesticides were eluted with high recoveries.

Since the target pesticides of this study ranged from relatively weak to strong acids, a silica gel column was used for further cleanup instead of a PSA or NH<sub>2</sub> anion exchange column, as in the Japanese official Multi-residue method II. As shown in Fig. 1b, the silica gel column cleanup procedure in the official Multi-residue method II involves sample loading with 2 mL of acetone/hexane/triethylamine (20:80:0.5 v/v/v), washing with an additional 10 mL of the same solvent mixture to remove low-polarity matrices, and then elution of pesticides with 10 mL of acetone/methanol (1:1 v/v). However, in this study, we omitted the washing procedure with acetone/hexane/triethylamine (20:80:0.5 v/v/v) because low-polarity pesticides elute together with low-polarity matrix components. Instead, we added the above-mentioned ODS column cleanup procedure to remove these matrices. Figure 3 shows the matrix effects of the representative pesticides in orange by eluting with various solvents from the silica gel column. Although eluting with acetone/methanol (1:1 v/v), the elution solvent used in the Japanese

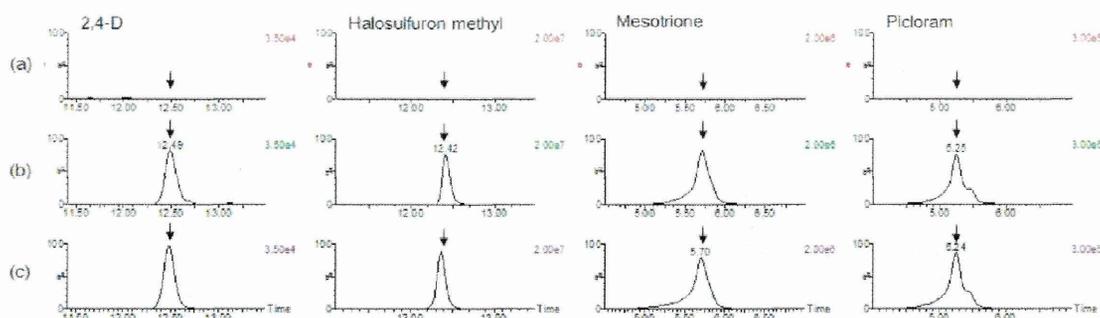


Fig. 2. Representative chromatograms of (a) extract of blank apple, (b) extract of apple spiked with  $0.01 \text{ mg kg}^{-1}$ , and (c)  $0.005 \mu\text{g mL}^{-1}$  standard solutions.

official Multi-residue method II achieved high recoveries of tested pesticides, large amounts of high-polarity components eluted along with the target pesticides, causing a significant matrix effect. In contrast, acetone alone eluted a small amount of matrix components from the silica gel column but yielded low recoveries of tested pesticides. Addition of formic acid to acetone improved pesticide recoveries but eluted matrix components and showed significant matrix effects. In contrast, addition of hexane, a lower-polarity solvent, to acetone/formic acid resulted in considerably decreased matrix effects for most of the target pesticides. Thus, by varying the proportion of acetone, hexane, and formic acid, the elution solvent was optimized. As a result, acetone/hexane/formic acid (100:100:1 v/v/v) has been shown to be an optimal elution solvent, giving minimal matrix effects with satisfactory recoveries except for mesotrione, nitenpyram, pyrasulfotole, and


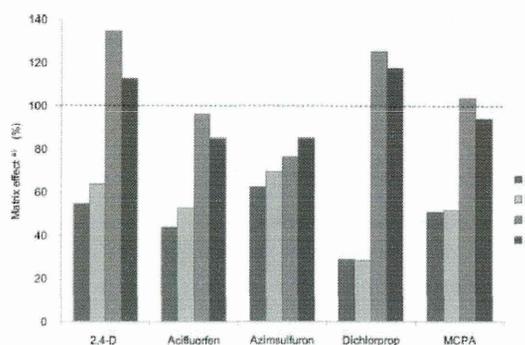


Fig. 3. Effect of elution solvents from silica gel column on matrix effects<sup>a</sup> of pesticide peak areas using an orange sample; a: methanol/acetone (1:1), 10 mL; b: acetone/formic acid (200:1), 10 mL; c: hexane/acetone/formic acid (50:150:1), 20 mL; d: hexane/acetone/formic acid (100:100:1), 20 mL. <sup>a</sup>Peak area ratio of a matrix-matched standard to standard in a pure solvent.

trinexapac. Although most of the pesticides were eluted with 10 mL of acetone/hexane/formic acid (100:100:1 v/v/v), foramsulfuron required 20 mL to elute. Thus, 20 mL of acetone/hexane/formic acid (100:100:1 v/v/v) was considered to give appropriate elution conditions for silica gel column cleanup.

Although the combination of ODS and silica gel column cleanup has been shown to be very effective for vegetables and fruits, the removal of pigments was not sufficient for samples containing high amounts of pigments, especially for spinach. Many methods for vegetables have been published using graphitized carbon sorbents for removing pigments.<sup>[8,9,22]</sup> However, for acidic pesticides, few methods using graphitized carbon have been reported, since acidic pesticides are prone to retention on the graphitized carbon sorbent. This can be explained by the anion exchange properties of positively charged oxonium groups, hydrophobic interactions between graphite carbon and the aromatic structure of pesticides, and hydrogen bonding between protonated functional groups of pesticides and carbonyl groups of the graphitized carbon or vice versa.<sup>[23]</sup> However, because graphitized carbon effectively removes pigments, we attempted to elute acidic pesticides from a graphitized carbon column by examining various solvents. Recoveries of representative pesticides from graphitized carbon column using different elution solvents are shown in Fig. 4. By using solvent standards, five pesticides (cyclanilide, flusulfamide, forchlorfenuron, pyrasulfotole, and teclotalam) out of 73 showed poor recoveries (<60%) by eluting with 22 mL of acetonitrile/toluene (3:1 v/v), which are the elution conditions used in the Japanese official Multi-residue method I. Recoveries of the tested pesticides changed slightly by varying the ratio of acetonitrile and toluene to (9:1 v/v) or (1:1 v/v). Moreover, the recoveries of most of the pesticides were affected only minimally by the addition of formic acid to acetonitrile/toluene (3:1 v/v). In contrast, addition of water to acetonitrile/toluene (3:1 v/v) gave significantly better recoveries. Thus, the ratio of water was optimized by comparing the recovery of

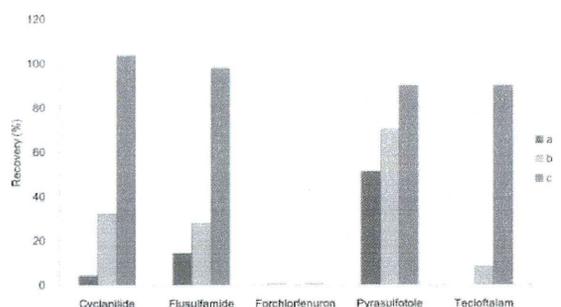


**Fig. 4.** Effect of elution solvents on the recoveries of pesticides from graphitized carbon column: a: acetonitrile/toluene (3:1), 20 mL; b: acetonitrile/toluene (1:1), 20 mL; c: acetonitrile/toluene/water (30:10:1), 20 mL.

pesticides from the graphitized carbon column with an acetonitrile/toluene/water eluting solvent of volume compositions of 30:10:0.2 v/v/v, 30:10:1 v/v/v, and 30:10:2 v/v/v. Considerably higher recoveries were obtained by eluting with acetonitrile/toluene/water 30:10:1 v/v/v compared with 30:10:0.2 v/v/v, while no great difference in

recoveries was observed between eluting with acetonitrile/toluene/water 30:10:1 v/v/v and 30:10:2 v/v/v. Most of the tested pesticides showed sufficient recoveries with 15 mL of acetonitrile/toluene/water (30:10:1 v/v/v), but flusulfamide and teclotalam required 25 mL to elute from the column. Only forchlorfenuron (0%) and trinexapac (66%) showed poor recoveries by eluting with 25 mL of acetonitrile/toluene/water (30:10:1 v/v/v). Therefore, 25 mL of acetonitrile/toluene/water (30:10:1 v/v/v) was considered to be the optimal elution solvent composition for use with the graphitized carbon column and can achieve high recoveries for most of the tested pesticides. The additional graphitized carbon column cleanup procedure effectively removed pigments, yielding a colorless solution from the spinach extract. The overall scheme of the above-described analytical procedure is shown in Fig. 5.

#### Recovery test

The recovery tests were performed for five times for each sample of cabbage, spinach, potato, eggplant, orange, and

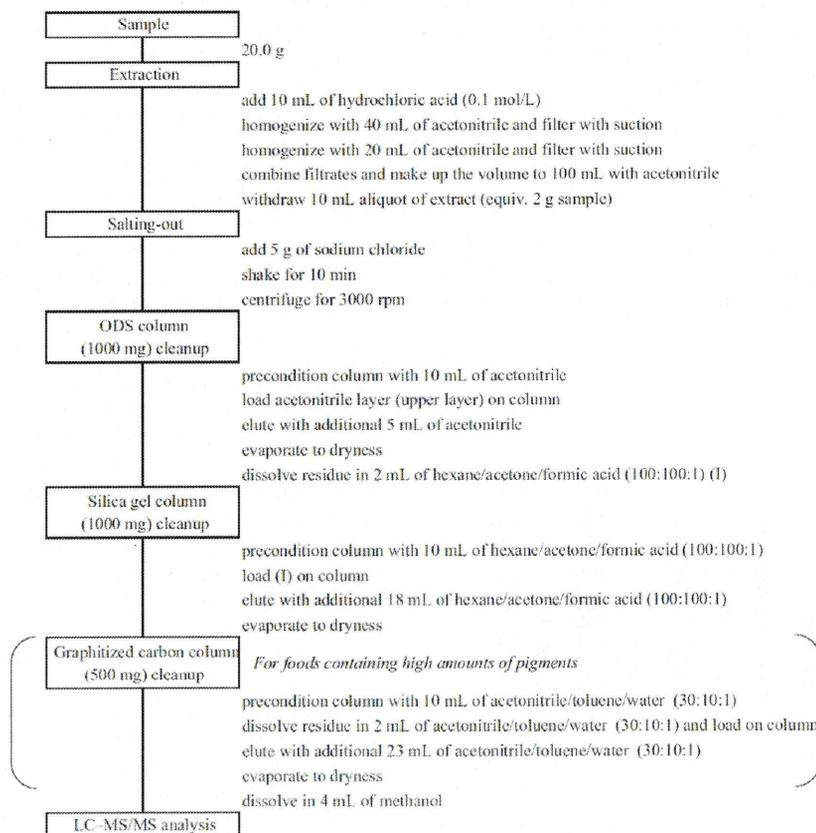

**Fig. 5.** Flow chart showing the sequence of steps in the developed multi-residue method.

Table 1. Recoveries of pesticides from fortified vegetables and fruits.

|                        | Cabbage      |         | Spinach      |         | Potato       |         | Eggplant     |         | Orange       |         | Apple        |         |
|------------------------|--------------|---------|--------------|---------|--------------|---------|--------------|---------|--------------|---------|--------------|---------|
|                        | Recovery (%) | RSD (%) |
| 2,4-D                  | 84           | 6       | 76           | 6       | 88           | 7       | 85           | 10      | 87           | 7       | 87           | 11      |
| 2,4-DB                 | 95           | 5       | 75           | 23      | 77           | 8       | 73           | 12      | 94           | 17      | 93           | 14      |
| 4-CPA                  | 84           | 10      | 70           | 13      | 91           | 16      | 82           | 12      | 60           | 14      | 94           | 14      |
| Acifluorfen            | 95           | 4       | 72           | 11      | 82           | 12      | 73           | 12      | 83           | 10      | 88           | 11      |
| Azimsulfuron           | 84           | 2       | 85           | 3       | 87           | 1       | 80           | 3       | 81           | 4       | 93           | 2       |
| Bensulfuron methyl     | 93           | 1       | 88           | 3       | 97           | 1       | 93           | 3       | 99           | 3       | 89           | 2       |
| Bentazone              | 85           | 3       | 77           | 6       | 86           | 5       | 86           | 5       | 21           | 22      | 85           | 7       |
| Bispyribac sodium      | 83           | 2       | 82           | 2       | 94           | 1       | 94           | 3       | 91           | 2       | 82           | 2       |
| Bromoxynil             | 94           | 5       | 80           | 13      | 76           | 7       | 75           | 4       | 83           | 10      | 89           | 9       |
| Bupirimate             | 87           | 1       | 84           | 2       | 87           | 2       | 84           | 3       | 81           | 2       | 77           | 3       |
| Chlorimuron ethyl      | 88           | 2       | 84           | 4       | 95           | 1       | 83           | 2       | 102          | 3       | 85           | 3       |
| Chlorsulfuron          | 86           | 1       | 83           | 4       | 88           | 2       | 88           | 3       | 96           | 1       | 88           | 2       |
| Cinosulfuron           | 82           | 2       | 89           | 3       | 92           | 1       | 90           | 4       | 105          | 2       | 94           | 1       |
| Clodinafop acid        | 80           | 1       | 86           | 2       | 72           | 3       | 72           | 4       | 91           | 2       | 79           | 3       |
| Cloprop                | 91           | 5       | 74           | 23      | 87           | 12      | 78           | 16      | 101          | 16      | 91           | 18      |
| Cloransulam methyl     | 92           | 2       | 90           | 5       | 101          | 2       | 92           | 3       | 114          | 3       | 100          | 2       |
| Cyclanilide            | 99           | 3       | 75           | 13      | 72           | 6       | 84           | 6       | 87           | 9       | 89           | 7       |
| Cyclosulfamuron        | 89           | 2       | 87           | 2       | 91           | 2       | 83           | 3       | 79           | 3       | 87           | 3       |
| Dichlorprop            | 83           | 9       | 62           | 16      | 100          | 10      | 102          | 5       | 86           | 8       | 88           | 10      |
| Diclosulam             | 88           | 1       | 88           | 3       | 94           | 1       | 90           | 2       | 81           | 2       | 85           | 2       |
| Ethametsulfuron methyl | 87           | 3       | 84           | 3       | 91           | 2       | 82           | 3       | 92           | 3       | 89           | 2       |
| Ethoxysulfuron         | 84           | 2       | 83           | 3       | 89           | 1       | 81           | 3       | 83           | 2       | 81           | 3       |
| Fenoprop               | 103          | 6       | 81           | 6       | 75           | 3       | 96           | 8       | 99           | 8       | 88           | 9       |
| Flazasulfuron          | 83           | 2       | 79           | 2       | 75           | 4       | 71           | 2       | 81           | 3       | 76           | 2       |
| Florasulam             | 90           | 0.4     | 86           | 4       | 89           | 1       | 81           | 1       | 114          | 1       | 83           | 1       |
| Fluazifop              | 92           | 3       | 72           | 2       | 74           | 3       | 85           | 2       | 97           | 1       | 85           | 2       |
| Flucarbazone sodium    | 91           | 1       | 87           | 4       | 100          | 1       | 86           | 2       | 96           | 3       | 92           | 3       |
| Flucetosulfuron        | 95           | 1       | 112          | 4       | 116          | 1       | 86           | 3       | 88           | 1       | 101          | 2       |
| Flumetsulam            | 92           | 1       | 88           | 4       | 92           | 1       | 73           | 3       | 86           | 2       | 87           | 2       |
| Flusulfamide           | 75           | 9       | 85           | 7       | 84           | 3       | 76           | 4       | 62           | 4       | 81           | 8       |
| Fomesafen              | 89           | 3       | 78           | 2       | 75           | 2       | 79           | 9       | 74           | 3       | 85           | 3       |
| Foramsulfuron          | 83           | 2       | 83           | 4       | 89           | 0.5     | 75           | 4       | 93           | 3       | 78           | 4       |
| Forchlorfenuron        | 89           | 2       | 14           | 17      | 84           | 3       | 29           | 3       | 74           | 1       | 71           | 2       |
| Gibberellie acid       | 78           | 2       | 76           | 4       | 81           | 5       | 81           | 4       | 66           | 3       | 83           | 6       |
| Halosulfuron methyl    | 74           | 2       | 83           | 4       | 86           | 1       | 64           | 1       | 77           | 3       | 72           | 4       |
| Haloxifop              | 87           | 2       | 83           | 2       | 91           | 5       | 82           | 3       | 73           | 2       | 87           | 3       |
| Imazamethabenz methyl  | 84           | 1       | 84           | 1       | 76           | 2       | 79           | 3       | 70           | 3       | 87           | 1       |
| Imazaquin              | 85           | 2       | 82           | 3       | 82           | 3       | 86           | 2       | 70           | 3       | 87           | 3       |
| Imazosulfuron          | 80           | 2       | 74           | 4       | 72           | 3       | 71           | 5       | 94           | 4       | 81           | 2       |
| Iodosulfuron methyl    | 89           | 1       | 91           | 3       | 106          | 2       | 91           | 3       | 93           | 1       | 77           | 4       |
| Ioxynil                | 85           | 4       | 73           | 3       | 80           | 3       | 82           | 8       | 50           | 17      | 91           | 4       |
| MCPA                   | 89           | 4       | 86           | 18      | 95           | 8       | 89           | 10      | 102          | 8       | 77           | 5       |
| Mecoprop               | 77           | 5       | 84           | 21      | 96           | 11      | 101          | 11      | 81           | 15      | 89           | 7       |
| Mesosulfuron methyl    | 84           | 2       | 93           | 4       | 99           | 1       | 85           | 4       | 119          | 2       | 92           | 2       |
| Mesotrione             | 78           | 2       | 76           | 5       | 74           | 3       | 80           | 2       | 103          | 2       | 92           | 2       |
| Metosulam              | 87           | 1       | 93           | 5       | 97           | 2       | 77           | 1       | 107          | 3       | 91           | 2       |
| Metsulfuron methyl     | 81           | 2       | 88           | 4       | 99           | 1       | 99           | 3       | 119          | 3       | 89           | 3       |
| Naptalam               | 72           | 1       | 71           | 4       | 73           | 2       | 73           | 1       | 75           | 6       | 72           | 3       |
| Nitenpyram             | 0            | —       | 0            | —       | 0            | —       | 0            | —       | 0            | —       | 0            | —       |
| Oryzalin               | 89           | 3       | 84           | 2       | 75           | 6       | 71           | 3       | 70           | 2       | 85           | 4       |
| Penoxsulam             | 92           | 1       | 94           | 3       | 101          | 2       | 92           | 1       | 105          | 2       | 97           | 2       |
| Picloram               | 78           | 3       | 58           | 2       | 60           | 3       | 82           | 6       | 79           | 11      | 83           | 2       |

(Continued on next page)

Table 1. Recoveries of pesticides from fortified vegetables and fruits. (Continued)

|                       | Cabbage      |         | Spinach      |         | Potato       |         | Eggplant     |         | Orange       |         | Apple        |         |
|-----------------------|--------------|---------|--------------|---------|--------------|---------|--------------|---------|--------------|---------|--------------|---------|
|                       | Recovery (%) | RSD (%) |
| Primisulfuron methyl  | 72           | 8       | 88           | 20      | 109          | 10      | 82           | 10      | 76           | 8       | 81           | 8       |
| Propoxycarbazone      | 92           | 1       | 88           | 2       | 91           | 1       | 83           | 2       | 114          | 0.4     | 97           | 2       |
| Propyrisulfuron       | 75           | 2       | 85           | 3       | 94           | 1       | 88           | 3       | 80           | 2       | 83           | 3       |
| Prosulfuron           | 83           | 6       | 86           | 2       | 97           | 2       | 79           | 3       | 96           | 4       | 76           | 4       |
| Pyrasulfotole         | 0            | —       | 0            | —       | 0            | —       | 0            | —       | 0            | —       | 0            | —       |
| Pyrazosulfuron ethyl  | 83           | 1       | 84           | 4       | 88           | 1       | 80           | 2       | 92           | 3       | 74           | 5       |
| Pyriothiobac sodium   | 89           | 3       | 77           | 2       | 75           | 2       | 85           | 3       | 87           | 1       | 73           | 4       |
| Quinclorac            | 79           | 3       | 70           | 3       | 70           | 3       | 75           | 3       | 73           | 5       | 84           | 1       |
| Rimsulfuron           | 80           | 2       | 84           | 3       | 83           | 4       | 82           | 4       | 83           | 4       | 85           | 2       |
| Sulfosulfuron         | 90           | 2       | 86           | 3       | 75           | 2       | 82           | 3       | 96           | 3       | 92           | 3       |
| Tecloftalam           | 36           | 13      | 50           | 16      | 47           | 5       | 43           | 8       | 41           | 10      | 76           | 6       |
| Thifensulfuron methyl | 78           | 1       | 90           | 3       | 99           | 1       | 99           | 3       | 110          | 2       | 96           | 1       |
| Tralkoxydim           | 79           | 4       | 83           | 3       | 70           | 4       | 85           | 4       | 79           | 4       | 81           | 3       |
| Triasulfuron          | 94           | 1       | 85           | 4       | 95           | 1       | 90           | 3       | 96           | 2       | 91           | 2       |
| Tribenuron methyl     | 77           | 6       | 86           | 3       | 78           | 7       | 43           | 7       | 58           | 8       | 47           | 8       |
| Triclopyr             | 90           | 7       | 104          | 10      | 78           | 7       | 83           | 8       | 81           | 7       | 73           | 12      |
| Trifloxysulfuron      | 84           | 3       | 87           | 3       | 88           | 5       | 78           | 4       | 90           | 5       | 77           | 5       |
| Triflusulfuron methyl | 104          | 3       | 87           | 3       | 97           | 2       | 96           | 4       | 96           | 3       | 91           | 4       |
| Trinexapac            | 95           | 6       | 71           | 5       | 73           | 1       | 75           | 3       | 77           | 4       | 92           | 4       |
| Trinexapac ethyl      | 84           | 2       | 80           | 2       | 81           | 4       | 73           | 1       | 86           | 5       | 86           | 1       |
| Warfarin <sup>a</sup> | 87           | 2       | 78           | 5       | 77           | 1       | 83           | 6       | 76           | 1       | 88           | 3       |

<sup>a</sup>Fortified at 0.001 mg kg<sup>-1</sup>.

apple fortified at 0.01 mg kg<sup>-1</sup>, except for warfarin at 0.001 mg kg<sup>-1</sup>. A spike level of 0.01 mg kg<sup>-1</sup> was chosen because, in Japan, a uniform limit of 0.01 mg kg<sup>-1</sup> is applied to all pesticides for which MRLs are not established. In addition, MRLs of warfarin have been set at 0.001 mg kg<sup>-1</sup> in various foods. Although matrix-matched calibration is often used to compensate for matrix effects, it is not necessarily accurate in real-world analyses because matrix effects vary with different blank samples. Thus, external calibration using solvent standards, instead of matrix-matched calibration, was used for quantifications. The results of the recovery tests are shown in Table 1. Out of the 73 tested pesticides, 70 for cabbage, 67 for spinach, 69 for potato, 67 for eggplant, 64 for orange, and 70 for apple were within 70–120%, with relative standard deviations below 25%, and were within the acceptable range of the Japanese guideline value.<sup>[24]</sup> Nitenpyram and pyrasulfotole were not recovered in all the samples tested. A possible reason for low recoveries of nitenpyram and pyrasulfotole would be their low recovery from the silica gel column, while poor recoveries of forchlorfenuron for the spinach sample can be explained by its insufficient elution from the graphitized carbon column. Trinexapac and mesotrione achieved high recoveries, in spite of poor recoveries in pure solvents from the silica gel column, indicating that coexisting matrix components facilitated the elution of these pesticides. The matrix effects were

evaluated by comparing the peak areas of matrix-matched standards with those of the standards in solvent (Table A2). Although some pesticides, such as picloram, tecloftalam, and tribenuron methyl, gave poor recoveries by ion suppression, the test solutions obtained by the developed method were much cleaner compared with those obtained by the Japanese official Multi-residue method II. A satisfactory linearity for all the tested pesticides was obtained with a correlation coefficient of >0.995. No significant interfering peak was observed in all blank samples, which indicated high selectivity of the developed method.

## Conclusions

In this study, we developed a sensitive and reliable LC-MS/MS method for the determination of a wide range of acidic pesticides in vegetables and fruits by modifying the Japanese official Multi-residue method II. The combination of ODS, silica gel, and graphitized carbon column cleanups effectively removed co-extracted matrix components and was suitable for subsequent LC-MS/MS analysis. Recoveries and precisions of most of the 73 tested pesticides were within the acceptable range of the Japanese guideline. Overall, the results suggest that the newly developed method is suitable for routine monitoring of acidic pesticide residues in vegetables and fruits.

## Funding

This work was supported by a Health and Labour Sciences Research Grant from the Ministry of Health, Labour, and Welfare in Japan.

## References

- [1] Mostafalou, S.; Abdollahi, M. Pesticides and human chronic diseases: evidences, mechanisms, and perspectives. *Toxicol. Appl. Pharmacol.* **2013**, *268*, 157–177.
- [2] Köhler, H.R.; Triebkorn, R. Wildlife ecotoxicology of pesticides: can we track effects to the population level and beyond? *Science* **2013**, *341*, 759–765.
- [3] Anastassiades, M.; Lehotay, S.J.; Stajnbaher, D.; Schenck, F.J. Fast and easy multi-residue method employing acetonitrile extraction/partitioning and “dispersive solid-phase extraction” for the determination of pesticide residues in produce. *J. AOAC Int.* **2003**, *86*, 412–431.
- [4] Lehotay, S.J.; Mastovská, K.; Lightfield, A.R. Use of buffering and other means to improve results of problematic pesticides in a fast and easy method for residue analysis of fruits and vegetables. *J. AOAC Int.* **2005**, *88*, 615–629.
- [5] Payá, P.; Anastassiades, M.; Mack, D.; Sigalova, I.; Tasdelen, B.; Oliva, J.; Barba, A. Analysis of pesticide residues using the quick easy cheap effective rugged and safe (QuEChERS) pesticide multi-residue method in combination with gas and liquid chromatography and tandem mass spectrometric detection. *Anal. Bioanal. Chem.* **2007**, *389*, 1697–1714.
- [6] Lehotay, S.J.; Mastovská, K.; Yun, S.J. Evaluation of two fast and easy methods for pesticide residue analysis in fatty food matrices. *J. AOAC Int.* **2005**, *88*, 630–638.
- [7] Koesukwiwat, U.; Lehotay, S.J.; Mastovska, K.; Dorweiler, K.J.; Leepipatpiboon, N. Extension of the QuEChERS method for pesticide residues in cereals to flaxseeds, peanuts, and doughs. *J. Agric. Food Chem.* **2010**, *58*, 5950–5958.
- [8] Fillion, J.; Sauvé, F.; Selwyn, J. Multi-residue method for the determination of residues of 251 pesticides in fruits and vegetables by gas chromatography/mass spectrometry and liquid chromatography with fluorescence detection. *J. AOAC Int.* **2000**, *83*, 698–713.
- [9] Ministry of Health, Labour and Welfare. *Director Notice, Syokusan No. 1129002*. Ministry of Health, Labour and Welfare: Tokyo, Japan, 2005.
- [10] Santilio, A.; Stefanelli, P.; Dommarco, R. Fast determination of phenoxy acid herbicides in carrots and apples using liquid chromatography coupled triple quadrupole mass spectrometry. *J. Environ. Sci. Health B* **2009**, *44*, 584–590.
- [11] Schaner, A.; Konecny, J.; Luckey, L.; Hickes, H. Determination of chlorinated acid herbicides in vegetation and soil by liquid chromatography/electrospray–tandem mass spectrometry. *J. AOAC Int.* **2007**, *90*, 1402–1410.
- [12] Koesukwiwat, U.; Sanguankaew, K.; Leepipatpiboon, N. Rapid determination of phenoxy acid residues in rice by modified QuEChERS extraction and liquid chromatography–tandem mass spectrometry. *Anal. Chim. Acta* **2008**, *626*, 10–20.
- [13] Charlton, A.J.A.; Stuckey, V.; Sykes, M.D. Determination of phenoxy acid herbicides MCPA, mecoprop and 2,4-D in kidney tissue using liquid chromatography with electrospray tandem mass spectrometry. *Bull. Environ. Contam. Toxicol.* **2009**, *82*, 711–715.
- [14] Kang, S.; Chang, N.; Zhao, Y.; Pan, C. Development of a method for the simultaneous determination of six sulfonylurea herbicides in wheat, rice, and corn by liquid chromatography–tandem mass spectrometry. *J. Agric. Food Chem.* **2011**, *59*, 9776–9781.
- [15] She, Y.X.; Cao, W.Q.; Shi, X.M.; Lv, X.L.; Liu, J.J.; Wang, R.Y.; Jin, F.; Wang, J.; Xiao, H. Class-specific molecularly imprinted polymers for the selective extraction and determination of sulfonylurea herbicides in maize samples by high-performance liquid chromatography–tandem mass spectrometry. *J. Chromatogr. B* **2010**, *878*, 2047–2053.
- [16] Peng, Y.; Xie, Y.; Luo, J.; Nie, L.; Chen, Y.; Chen, L.; Du, S.; Zhang, Z. Molecularly imprinted polymer layer-coated silica nanoparticles toward dispersive solid-phase extraction of trace sulfonylurea herbicides from soil and crop samples. *Anal. Chim. Acta* **2010**, *674*, 190–200.
- [17] Quesada-Molina, C.; del Olmo-Iruela, M.; García-Campana, A.M. Trace determination of sulfonylurea herbicides in water and grape samples by capillary zone electrophoresis using large volume sample stacking. *Anal. Bioanal. Chem.* **2010**, *397*, 2593–2601.
- [18] Seccia, S.; Albrizio, S.; Fidente, P.; Montesano, D. Development and validation of a solid-phase extraction method coupled to high-performance liquid chromatography with ultraviolet–diode array detection for the determination of sulfonylurea herbicide residues in bovine milk samples. *J. Chromatogr. A* **2011**, *1218*, 1253–1259.
- [19] Pareja, L.; Ccsio, V.; Heinzen, H.; Fernández-Alba, A.R. Evaluation of various QuEChERS-based methods for the analysis of herbicides and other commonly used pesticides in polished rice by LC–MS/MS. *Talanta* **2011**, *83*, 1613–1622.
- [20] Akiyama, Y.; Matsuoka, T.; Mitsuhashi, T. Multi-residue screening method of acidic pesticides in agricultural products by liquid chromatography/time of flight mass spectrometry. *J. Pestic. Sci.* **2009**, *34*, 265–272.
- [21] Okihashi, M.; Kitagawa, Y.; Akutsu, K.; Obana, H.; Tanaka, Y. Rapid method for the determination of 180 pesticide residues in foods by gas chromatography/mass spectrometry and flame photometric detection. *J. Pestic. Sci.* **2005**, *30*, 368–377.
- [22] Zheng, H.B.; Zhao, Q.; Mo J.Z.; Huang, Y.Q.; Luo, Y.B.; Yu, Q.W.; Feng, Y.Q. Quick, easy, cheap, effective, rugged and safe method with magnetic graphitized carbon black and primary secondary amine as adsorbent and its application in pesticide residue analysis. *J. Chromatogr. A* **2013**, *1300*, 127–133.
- [23] Laganà, A.; Bacaloni, A.; De Leva, I.; Faberi, A.; Fago, G.; Marino, A. Occurrence and determination of herbicides and their major transformation products in environmental waters. *Anal. Chim. Acta* **2002**, *462*, 187–198.
- [24] Ministry of Health, Labour and Welfare. *Director Notice, Syokusan No. 1224-1*. Ministry of Health, Labour and Welfare: Tokyo, Japan, 2010.

## Appendix

Table A1. LC-MS/MS parameters for the tested pesticides.

|                        | Ionization mode | Retention time (min) | Quantification      |                   |                  |                       | Confirmation        |                   |                  |                       |
|------------------------|-----------------|----------------------|---------------------|-------------------|------------------|-----------------------|---------------------|-------------------|------------------|-----------------------|
|                        |                 |                      | Precursor ion (m/z) | Product ion (m/z) | Cone voltage (V) | Collision energy (eV) | Precursor ion (m/z) | Product ion (m/z) | Cone voltage (V) | Collision energy (eV) |
| 2,4-D                  | ESI(-)          | 12.5                 | 219.0               | 160.9             | 30               | 10                    | 219.0               | 125.1             | 30               | 25                    |
| 2,4-DB                 | ESI(-)          | 15.3                 | 160.8               | 124.9             | 40               | 15                    | 247.1               | 161.0             | 60               | 15                    |
| 4-CPA                  | ESI(-)          | 9.9                  | 185.0               | 127.0             | 40               | 15                    | 187.0               | 127.0             | 40               | 15                    |
| Acifluorfen            | ESI(-)          | 16.1                 | 360.1               | 316.0             | 10               | 10                    | 360.1               | 195.0             | 10               | 25                    |
| Azimsulfuron           | ESI(+)          | 9.8                  | 425.3               | 182.2             | 50               | 20                    | 425.3               | 156.2             | 50               | 35                    |
| Bensulfuron methyl     | ESI(+)          | 15.0                 | 411.3               | 149.2             | 50               | 20                    | 411.3               | 182.2             | 50               | 20                    |
| Bentazone              | ESI(-)          | 8.5                  | 239.1               | 132.1             | 60               | 25                    | 239.1               | 197.0             | 60               | 20                    |
| Bispyribac sodium      | ESI(+)          | 15.1                 | 431.3               | 275.2             | 30               | 15                    | 431.3               | 413.2             | 30               | 15                    |
| Bromoxynil             | ESI(-)          | 10.8                 | 275.9               | 80.9              | 60               | 25                    | 275.9               | 78.9              | 60               | 25                    |
| Bupirimate             | ESI(+)          | 19.4                 | 317.3               | 166.2             | 50               | 25                    | 317.3               | 150.2             | 50               | 20                    |
| Chlorimuron ethyl      | ESI(+)          | 13.3                 | 415.2               | 186.1             | 40               | 20                    | 415.2               | 121.1             | 40               | 40                    |
| Chlorsulfuron          | ESI(+)          | 9.9                  | 358.2               | 141.2             | 40               | 15                    | 358.2               | 167.2             | 40               | 15                    |
| Cinosulfuron           | ESI(+)          | 9.9                  | 414.3               | 183.2             | 50               | 15                    | 414.3               | 157.2             | 50               | 25                    |
| Clodinafop acid        | ESI(+)          | 14.6                 | 312.2               | 266.2             | 50               | 15                    | 312.2               | 91.1              | 50               | 25                    |
| Cloprop                | ESI(-)          | 11.3                 | 199.0               | 127.0             | 10               | 15                    | 199.0               | 71.0              | 10               | 10                    |
| Cloransulam methyl     | ESI(+)          | 13.3                 | 430.2               | 398.1             | 40               | 15                    | 430.2               | 370.1             | 40               | 20                    |
| Cyclanilide            | ESI(-)          | 14.8                 | 272.0               | 160.0             | 40               | 20                    | 272.0               | 228.0             | 40               | 10                    |
| Cyclosulfamuron        | ESI(+)          | 16.0                 | 422.3               | 261.2             | 40               | 15                    | 422.3               | 218.2             | 40               | 25                    |
| Dichlorprop            | ESI(-)          | 14.0                 | 233.0               | 161.0             | 10               | 15                    | 233.0               | 125.0             | 10               | 25                    |
| Diclosulam             | ESI(+)          | 13.6                 | 406.1               | 161.1             | 40               | 25                    | 406.1               | 378.1             | 40               | 15                    |
| Ethametsulfuron methyl | ESI(+)          | 11.7                 | 411.3               | 196.2             | 40               | 15                    | 411.3               | 168.2             | 40               | 30                    |
| Ethoxysulfuron         | ESI(+)          | 13.6                 | 399.2               | 261.2             | 50               | 15                    | 399.2               | 218.2             | 50               | 25                    |
| Fenoprop               | ESI(-)          | 15.7                 | 266.8               | 194.7             | 10               | 15                    | 266.8               | 158.7             | 10               | 30                    |
| Flazasulfuron          | ESI(+)          | 10.7                 | 408.2               | 182.2             | 40               | 20                    | 408.2               | 139.2             | 40               | 40                    |
| Florasulam             | ESI(+)          | 9.8                  | 360.2               | 129.2             | 40               | 20                    | 360.2               | 192.1             | 40               | 15                    |
| Fluazifop              | ESI(+)          | 14.6                 | 328.2               | 282.2             | 50               | 20                    | 328.2               | 254.2             | 50               | 25                    |
| Flucarbazone sodium    | ESI(+)          | 10.4                 | 397.2               | 130.2             | 30               | 10                    | 397.2               | 115.2             | 30               | 40                    |
| Flucetosulfuron        | ESI(+)          | 11.3                 | 488.3               | 156.2             | 40               | 20                    | 488.3               | 273.2             | 40               | 25                    |
| Flumetsulam            | ESI(+)          | 7.5                  | 326.2               | 129.2             | 50               | 25                    | 326.2               | 262.2             | 50               | 20                    |
| Flusulfamide           | ESI(-)          | 18.2                 | 413.0               | 171.0             | 80               | 35                    | 413.0               | 349.0             | 80               | 25                    |
| Fomesafen              | ESI(+)          | 16.1                 | 456.2               | 344.1             | 30               | 15                    | 456.2               | 223.2             | 30               | 30                    |
| Foramsulfuron          | ESI(+)          | 10.2                 | 453.3               | 182.2             | 30               | 20                    | 453.3               | 272.2             | 30               | 10                    |
| Forchlorfenuron        | ESI(+)          | 16.4                 | 248.2               | 129.1             | 40               | 20                    | 248.2               | 93.1              | 40               | 30                    |
| Gibberellic acid       | ESI(+)          | 8.3                  | 364.3               | 239.3             | 30               | 15                    | 364.3               | 311.3             | 30               | 15                    |
| Halosulfuron methyl    | ESI(+)          | 12.4                 | 435.2               | 182.2             | 40               | 20                    | 435.2               | 139.2             | 40               | 40                    |
| Haloxypop              | ESI(+)          | 16.7                 | 362.2               | 288.2             | 40               | 25                    | 362.2               | 316.2             | 40               | 15                    |
| Imazamethabenz methyl  | ESI(+)          | 14.7                 | 289.3               | 144.1             | 50               | 35                    | 289.3               | 161.2             | 50               | 25                    |
| Imazaquin              | ESI(+)          | 9.7                  | 312.3               | 267.2             | 40               | 20                    | 312.3               | 199.2             | 40               | 25                    |
| Imazosulfuron          | ESI(+)          | 11.0                 | 413.2               | 153.2             | 40               | 10                    | 413.2               | 156.2             | 40               | 20                    |
| Iodosulfuron methyl    | ESI(+)          | 12.2                 | 508.1               | 167.2             | 40               | 20                    | 508.1               | 141.2             | 40               | 25                    |
| Joxynil                | ESI(-)          | 12.8                 | 369.9               | 127.0             | 40               | 30                    | 369.9               | 215.0             | 40               | 30                    |
| MCPA                   | ESI(-)          | 12.5                 | 199.1               | 141.0             | 30               | 15                    | 199.1               | 105.0             | 30               | 25                    |
| Mecoprop               | ESI(-)          | 13.9                 | 213.1               | 141.0             | 30               | 20                    | 213.1               | 71.1              | 30               | 10                    |
| Mesosulfuron methyl    | ESI(+)          | 11.4                 | 504.2               | 182.2             | 50               | 20                    | 504.2               | 162.2             | 50               | 40                    |
| Mesotrione             | ESI(+)          | 5.7                  | 340.2               | 228.1             | 40               | 15                    | 340.2               | 104.1             | 40               | 30                    |
| Metosulam              | ESI(+)          | 12.9                 | 418.2               | 175.2             | 50               | 25                    | 418.2               | 140.2             | 50               | 40                    |
| Metsulfuron methyl     | ESI(+)          | 8.8                  | 382.2               | 167.2             | 40               | 15                    | 382.2               | 199.2             | 40               | 20                    |

(Continued on next page)

Table A1. LC-MS/MS parameters for the tested pesticides. (Continued)

|                       | Ionization mode | Retention time (min) | Quantification      |                   |                  |                       | Confirmation        |                   |                  |                       |
|-----------------------|-----------------|----------------------|---------------------|-------------------|------------------|-----------------------|---------------------|-------------------|------------------|-----------------------|
|                       |                 |                      | Precursor ion (m/z) | Product ion (m/z) | Cone voltage (V) | Collision energy (eV) | Precursor ion (m/z) | Product ion (m/z) | Cone voltage (V) | Collision energy (eV) |
| Naptalam              | ESI(+)          | 12.1                 | 292.3               | 144.1             | 30               | 10                    | 292.3               | 149.1             | 30               | 20                    |
| Nitenpyram            | ESI(+)          | 7.0                  | 271.2               | 126.1             | 30               | 25                    | 271.2               | 130.2             | 30               | 10                    |
| Oryzalin              | ESI(+)          | 18.5                 | 347.2               | 288.2             | 40               | 15                    | 347.2               | 305.2             | 40               | 15                    |
| Penoxsulam            | ESI(+)          | 12.9                 | 484.2               | 195.2             | 60               | 25                    | 484.2               | 164.2             | 60               | 35                    |
| Picloram              | ESI(+)          | 5.2                  | 241.1               | 168.0             | 30               | 30                    | 241.1               | 195.0             | 30               | 20                    |
| Primisulfuron methyl  | ESI(-)          | 14.8                 | 467.1               | 226.1             | 10               | 15                    | 467.1               | 176.1             | 10               | 30                    |
| Propoxycarbazon       | ESI(+)          | 10.6                 | 399.1               | 198.9             | 20               | 10                    | 399.1               | 158.0             | 20               | 10                    |
| Propyrisulfuron       | ESI(+)          | 13.6                 | 456.2               | 261.2             | 40               | 15                    | 456.2               | 196.2             | 40               | 15                    |
| Prosulfuron           | ESI(+)          | 13.9                 | 420.2               | 141.2             | 50               | 20                    | 420.2               | 167.2             | 50               | 20                    |
| Pyrasulfotole         | ESI(+)          | 10.7                 | 363.2               | 251.1             | 50               | 25                    | 363.2               | 220.1             | 50               | 40                    |
| Pyrazosulfuron ethyl  | ESI(+)          | 12.1                 | 415.2               | 182.2             | 40               | 25                    | 415.2               | 139.2             | 40               | 40                    |
| Pyriithiobac sodium   | ESI(+)          | 13.6                 | 327.2               | 309.1             | 40               | 15                    | 327.2               | 139.2             | 40               | 25                    |
| Quinclorac            | ESI(+)          | 8.4                  | 242.1               | 161.1             | 30               | 35                    | 242.1               | 224.1             | 30               | 15                    |
| Rimsulfuron           | ESI(+)          | 9.1                  | 432.2               | 182.2             | 40               | 25                    | 432.2               | 325.2             | 40               | 15                    |
| Sulfosulfuron         | ESI(+)          | 10.7                 | 471.2               | 211.2             | 40               | 15                    | 471.2               | 261.2             | 40               | 15                    |
| Teclofitalam          | ESI(+)          | 17.9                 | 465.2               | 162.1             | 20               | 15                    | 466.9               | 162.1             | 20               | 15                    |
| Thifensulfuron methyl | ESI(+)          | 8.9                  | 388.2               | 167.2             | 40               | 15                    | 388.2               | 205.1             | 40               | 25                    |
| Tralkoxydim           | ESI(+)          | 15.4                 | 330.4               | 284.3             | 40               | 10                    | 330.4               | 138.2             | 40               | 20                    |
| Triasulfuron          | ESI(+)          | 10.9                 | 402.2               | 167.2             | 40               | 15                    | 402.2               | 141.2             | 40               | 20                    |
| Tribenuron methyl     | ESI(+)          | 11.1                 | 396.2               | 155.2             | 40               | 10                    | 396.2               | 181.2             | 40               | 20                    |
| Triclopyr             | ESI(+)          | 13.4                 | 256.1               | 146.0             | 30               | 25                    | 257.9               | 146.0             | 30               | 25                    |
| Trifloxysulfuron      | ESI(+)          | 12.2                 | 438.2               | 182.2             | 30               | 25                    | 438.2               | 139.2             | 30               | 40                    |
| Triflusulfuron methyl | ESI(+)          | 14.7                 | 493.3               | 264.2             | 40               | 20                    | 493.3               | 96.1              | 40               | 40                    |
| Trinexapac            | ESI(+)          | 3.1                  | 225.0               | 68.9              | 40               | 15                    | 225.0               | 164.9             | 40               | 15                    |
| Trinexapac ethyl      | ESI(+)          | 9.4                  | 253.0               | 207.0             | 20               | 10                    | 253.0               | 185.0             | 20               | 10                    |
| Warfarin              | ESI(+)          | 12.9                 | 309.3               | 163.2             | 30               | 15                    | 309.3               | 251.2             | 30               | 20                    |

Table A2. Matrix effects<sup>a</sup> of the tested pesticides.

|                    | Cabbage | Spinach | Potato | Eggplant | Orange | Apple |
|--------------------|---------|---------|--------|----------|--------|-------|
| 2,4-D              | 89      | 81      | 95     | 86       | 97     | 90    |
| 2,4-DB             | 100     | 100     | 78     | 76       | 111    | 83    |
| 4-CPA              | 85      | 73      | 82     | 76       | 56     | 88    |
| Acifluorfen        | 100     | 89      | 72     | 74       | 84     | 96    |
| Azimsulfuron       | 90      | 94      | 86     | 86       | 78     | 100   |
| Bensulfuron methyl | 100     | 97      | 103    | 102      | 105    | 95    |
| Bentazone          | 83      | 81      | 98     | 93       | 44     | 90    |
| Bispyribac sodium  | 97      | 97      | 104    | 107      | 99     | 92    |
| Bromoxynil         | 102     | 92      | 83     | 84       | 89     | 91    |
| Bupirimate         | 95      | 93      | 99     | 97       | 93     | 82    |
| Chlorimuron ethyl  | 96      | 94      | 100    | 91       | 105    | 87    |
| Chlorsulfuron      | 90      | 92      | 92     | 92       | 102    | 90    |
| Cinosulfuron       | 88      | 98      | 75     | 95       | 107    | 97    |
| Clodinafop acid    | 89      | 93      | 83     | 79       | 94     | 83    |
| Cloprop            | 97      | 72      | 112    | 105      | 106    | 88    |
| Cloransulam methyl | 100     | 99      | 106    | 102      | 116    | 101   |
| Cyclanilide        | 107     | 94      | 81     | 93       | 93     | 96    |
| Cyclosulfamuron    | 95      | 95      | 95     | 92       | 79     | 90    |
| Dichlorprop        | 97      | 68      | 99     | 102      | 99     | 100   |

(Continued on next page)

Table A2. Matrix effects<sup>a</sup> of the tested pesticides. (Continued)

|                        | <i>Cabbage</i> | <i>Spinach</i> | <i>Potato</i> | <i>Eggplant</i> | <i>Orange</i> | <i>Apple</i> |
|------------------------|----------------|----------------|---------------|-----------------|---------------|--------------|
| Diclosulam             | 95             | 96             | 101           | 100             | 87            | 86           |
| Ethametsulfuron methyl | 92             | 94             | 86            | 86              | 97            | 90           |
| Ethoxysulfuron         | 90             | 93             | 87            | 89              | 83            | 85           |
| Fenoprop               | 117            | 117            | 98            | 105             | 101           | 96           |
| Flazasulfuron          | 86             | 79             | 81            | 81              | 102           | 79           |
| Florasulam             | 95             | 95             | 94            | 93              | 112           | 92           |
| Fluazifop              | 92             | 73             | 80            | 94              | 98            | 91           |
| Flucarbazone sodium    | 99             | 100            | 105           | 94              | 100           | 95           |
| Flucetosulfuron        | 109            | 103            | 111           | 102             | 89            | 101          |
| Flumetsulam            | 99             | 97             | 100           | 96              | 92            | 97           |
| Flusulfamide           | 81             | 95             | 99            | 88              | 69            | 81           |
| Fomesafen              | 90             | 90             | 86            | 76              | 83            | 89           |
| Foramsulfuron          | 93             | 94             | 78            | 85              | 120           | 86           |
| Forchlorfenuron        | 95             | 61             | 98            | 45              | 80            | 77           |
| Gibberellic acid       | 84             | 80             | 97            | 92              | 78            | 95           |
| Halosulfuron methyl    | 83             | 91             | 90            | 72              | 94            | 78           |
| Haloxypop              | 92             | 92             | 98            | 96              | 76            | 89           |
| Imazamethabenz methyl  | 91             | 96             | 84            | 88              | 85            | 94           |
| Imazaquin              | 90             | 84             | 87            | 97              | 79            | 94           |
| Imazosulfuron          | 88             | 82             | 69            | 90              | 97            | 82           |
| Iodosulfuron methyl    | 96             | 103            | 110           | 96              | 105           | 79           |
| Ioxynil                | 90             | 78             | 88            | 90              | 68            | 94           |
| MCPA                   | 98             | 102            | 99            | 90              | 96            | 96           |
| Mecoprop               | 83             | 93             | 96            | 104             | 88            | 106          |
| Mesosulfuron methyl    | 89             | 102            | 100           | 92              | 118           | 92           |
| Mesotrione             | 94             | 95             | 95            | 95              | 111           | 96           |
| Metosulam              | 95             | 102            | 105           | 92              | 110           | 93           |
| Metsulfuron methyl     | 88             | 97             | 99            | 103             | 116           | 93           |
| Naptalam               | 92             | 91             | 96            | 99              | 97            | 88           |
| Nitenpyram             | 87             | 95             | 93            | 96              | 88            | 94           |
| Oryzalin               | 97             | 97             | 84            | 79              | 80            | 92           |
| Penoxsulam             | 100            | 105            | 104           | 103             | 108           | 98           |
| Picloram               | 86             | 64             | 69            | 98              | 93            | 84           |
| Primisulfuron methyl   | 94             | 105            | 120           | 75              | 94            | 86           |
| Propoxycarbazon        | 101            | 99             | 93            | 98              | 116           | 99           |
| Propyrisulfuron        | 83             | 92             | 86            | 95              | 85            | 83           |
| Prosulfuron            | 87             | 96             | 106           | 91              | 100           | 79           |
| Pyrasulfotole          | 92             | 89             | 82            | 59              | 86            | 89           |
| Pyrazosulfuron ethyl   | 92             | 94             | 82            | 87              | 93            | 71           |
| Pyriithiobac sodium    | 94             | 88             | 85            | 100             | 95            | 76           |
| Quinclorac             | 87             | 92             | 99            | 89              | 93            | 88           |
| Rimsulfuron            | 84             | 93             | 93            | 90              | 97            | 92           |
| Sulfosulfuron          | 98             | 97             | 73            | 93              | 106           | 95           |
| Tecloftalam            | 47             | 69             | 63            | 53              | 54            | 89           |
| Thifensulfuron methyl  | 86             | 97             | 101           | 101             | 116           | 93           |
| Tralkoxydim            | 90             | 90             | 78            | 98              | 86            | 84           |
| Triasulfuron           | 97             | 94             | 94            | 95              | 104           | 91           |
| Tribenuron methyl      | 92             | 75             | 95            | 52              | 67            | 54           |
| Triclopyr              | 82             | 84             | 79            | 84              | 79            | 79           |
| Trifloxysulfuron       | 91             | 97             | 79            | 88              | 90            | 85           |
| Triflusulfuron methyl  | 110            | 95             | 104           | 105             | 99            | 91           |
| Trinexapac             | 94             | 80             | 84            | 73              | 81            | 89           |
| Trinexapac ethyl       | 92             | 92             | 88            | 79              | 96            | 91           |
| Warfarin <sup>b</sup>  | 93             | 91             | 87            | 101             | 79            | 91           |

<sup>a</sup>Matrix effect (%) = peak area of matrix-matched standard/peak area of standard in solvent × 100.<sup>b</sup>Fortified at 0.001 mg/kg.

# Multiresidue determination of pesticides in tea by gas chromatography-tandem mass spectrometry

SHIZUKA SAITO-SHIDA, SATORU NEMOTO and REIKO TESHIMA

Division of Foods, National Institute of Health Sciences, Tokyo, Japan

An efficient and reliable GC-MS/MS method for the multiresidue determination of pesticides in tea was developed by modifying the Japanese official multiresidue method. Sample preparation was carefully optimized for the efficient removal of coextracted matrix components. The optimal sample preparation procedure involved swelling of the sample in water; extraction with acetonitrile; removal of water by salting-out; and sequential cleanup by ODS, graphitized carbon black/primary secondary amine (GCB/PSA) and silica gel cartridges prior to GC-MS/MS analysis. The recoveries of 162 pesticides from fortified (at 0.01 mg kg<sup>-1</sup>) green tea, oolong tea, black tea and matcha (powdered green tea) were mostly (95–98% of the tested pesticides) within the range of 70–120%, with relative standard deviations of <20%. Poor recovery of triazole pesticides was considered to be due to low recovery from the silica gel cartridges. The test solutions obtained by the modified method contained relatively small amounts of pigments, caffeine and other matrix components and were cleaner than those obtained by the original Japanese official multiresidue method. No interfering peaks were observed in the blank chromatograms, indicating the high selectivity of the modified method. The overall results suggest that the developed method is suitable for the quantitative analysis of GC-amenable pesticide residues in tea.

**Keywords:** GC-MS/MS, multiresidue method, pesticides, tea.

## Introduction

Tea (*Camellia sinensis*) is one of the most widely consumed beverages in the world. Depending on its level of fermentation during the manufacturing process, tea can be classified into several groups, including green tea (unfermented), oolong tea (partially fermented) and black tea (fully fermented).<sup>[1]</sup> In Japan, a finely ground powdered green tea, known as matcha, is also consumed. Although tea is known to have a number of health benefits such as anti-cancer properties,<sup>[2,3]</sup> the consumption of tea can also be a potential source of exposure to pesticides used for tea cultivation control, resulting in long-term health effects. In addition, the concentration of pesticide residues in tea is generally high compared to other foods. As tea leaves have a large surface area per mass and the mass of tea leaves is reduced during the drying process, tea has high pesticide concentration.<sup>[4]</sup>

For these reasons, there is a demand for reliable multiresidue methods for analyzing pesticides present in tea.

However, this is a challenging task, as tea leaves contain large amounts of complicated matrices such as pigments, polyphenols and caffeine, which contaminate LC-MS (/MS) and GC-MS (/MS) systems, thus interfering with analysis. Caffeine in particular is difficult to remove from tea extracts because it possesses similar chemical properties to a number of pesticides,<sup>[5]</sup> thus causing serious problems for GC-MS (/MS) analysis.

To date, there have been a number of reports on the multiresidue analysis of pesticide residues in tea using LC-MS (/MS) or GC-MS (/MS) techniques.<sup>[6–12]</sup> Recently, several modifications<sup>[13,14]</sup> of the QuEChERS (quick, easy, cheap, effective, rugged, safe) method, developed by Anastassiades et al.,<sup>[15]</sup> have been published. In Japan, an official “multiresidue method for agricultural chemicals” has been established (Fig. 1).<sup>[16]</sup> In addition, we previously reported the development of an LC-MS/MS method for the multiresidue analysis of pesticides in tea.<sup>[17]</sup> However, all these methods require frequent maintenance of the GC-MS/MS system, as they lack a cleanup step for the removal of caffeine from the system, resulting in a significant amount of caffeine in the later test solutions.

Thus, the aims of this study were as follows: (1) optimization of an efficient caffeine-removal cleanup step for the GC-MS/MS analysis of pesticide residues in tea; (2) development of a sensitive, selective and reliable analytical

Address correspondence to Shizuka Saito-Shida, Division of Foods, National Institute of Health Sciences, 1-18-1, Kamiyoga, Setagaya-ku, Tokyo 158-8501, Japan; E-mail: shizsaito@nih.go.jp

Received February 23, 2015.

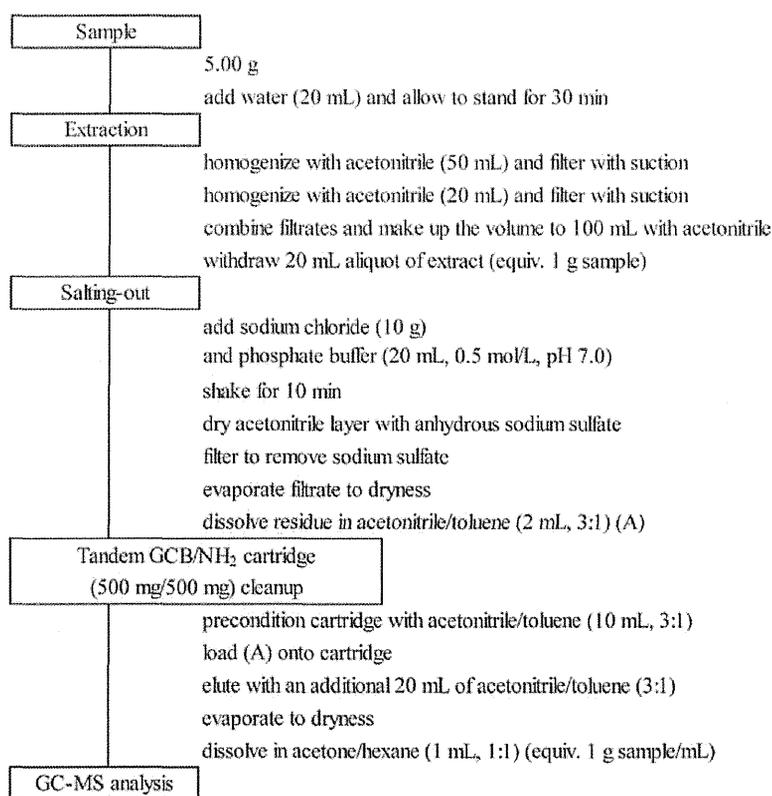


Fig. 1. Flow chart of sample preparation according to the Japanese official multiresidue method.

method for multiresidue determination of GC-amenable pesticides in tea by modifying the Japanese official multiresidue method; and (3) application of the proposed method to green tea, oolong tea, black tea and matcha.

## Experimental

### Reagents and materials

**Solvents and reagents.** Pesticide residue analysis-grade acetone, acetonitrile, hexane and toluene were obtained from Kanto Chemical (Tokyo, Japan). Diatomaceous earth (Celite, No. 545), pesticide residue analysis-grade sodium chloride, analytical-grade dipotassium hydrogen phosphate and analytical-grade potassium dihydrogen phosphate were purchased from Wako Pure Chemical Industries (Osaka, Japan). The water used to prepare the test solutions was purified using an NZJ-2DSYW distillation apparatus (Fujiwara Scientific, Tokyo, Japan).

**Analytical standards.** Pesticide standards were purchased from Dr. Ehrenstorfer (Augsburg, Germany), Riedel-de

Haën (Seelze, Germany), Sigma-Aldrich (St. Louis, MO, USA), Hayashi Pure Chemical (Osaka, Japan), Wako Pure Chemical Industries, and the Kanto Chemical Co. Individual stock standard solutions ( $1,000 \text{ mg L}^{-1}$ ) were prepared in hexane, or in a mixture of hexane and acetone, depending on the solubility of each pesticide. Working standard solutions were prepared by mixing the stock standard solutions, and diluting with acetone and hexane as required.

**Cartridge columns.** Octadecylsilyl silica gel (ODS) cartridges (Mega Bond Elut C18, 1,000 mg) were obtained from Agilent Technologies (Palo Alto, CA, USA). Tandem graphitized carbon black/primary secondary amine (GCB/PSA) cartridges (InertSep GC/PSA, 500 mg/500 mg) and silica gel cartridges (InertSep SI, 1,000 mg) were obtained from GL Sciences (Tokyo, Japan).

**Phosphate buffer.** Phosphate buffer ( $0.5 \text{ mol L}^{-1}$ , pH 7.0) was prepared by dissolving dipotassium hydrogen phosphate (52.7 g) and potassium dihydrogen phosphate (30.2 g) in water (500 mL). The pH was then adjusted to

pH 7.0 by the addition of either 1 mol L<sup>-1</sup> sodium hydroxide or 1 mol L<sup>-1</sup> hydrochloric acid. The final volume of the solution was made up to 1 L with water.

**Food samples.** Green tea leaves, oolong tea leaves, black tea leaves and matcha (powdered green tea), were purchased from a market in Tokyo (Japan). Green tea leaves, oolong tea leaves and black tea leaves were ground into small particles using a centrifugal mill and passed through a 425  $\mu$ m pore standard sieve.

#### Apparatus

**GC-MS/MS.** A Trace 1310 Gas Chromatograph (ThermoFisher Scientific, MA, USA) coupled to a TSQ 8000 Mass Spectrometer (ThermoFisher Scientific), equipped with a TriPlus RSH Autosampler (ThermoFisher Scientific) was used under the following operating conditions: column, DB-5ms column (30 m length, 0.25 mm id, 0.25  $\mu$ m film thickness; Agilent) equipped with a guard column (2 m length, 0.25 mm id; Agilent); column temperature, 50 °C (1 min), increased to 125 °C at 25 °C min<sup>-1</sup>, then increased to 300 °C at 10 °C min<sup>-1</sup> and held for 8.5 min; carrier gas, helium; flow rate, 1 mL min<sup>-1</sup>; injection volume, 2  $\mu$ L; ionization mode, electron impact (EI); ion source temperature, 260 °C; injection port temperature, 260 °C; transfer line temperature, 280 °C. The selected reaction monitoring (SRM) transitions and retention times are summarized in Table A1.

**Laboratory apparatus.** A homogenizer (Polytron PT 10–35 GT; Kinematica, Lucerne, Switzerland), centrifugal mill (Ultra Centrifugal Mill ZM 200; Retsch, Haan, Germany), rotary evaporator (N-1000/NVC-2100; Tokyo Rikakikai, Tokyo, Japan), electric shaker (SR-2w; Taitec, Saitama, Japan) and centrifuge (Centrifuge 8100; Kubota, Tokyo, Japan) were used in the studies as described.

#### Sample preparation

The desired sample (5.00 g) of tea was weighed in a 350 mL glass tube and was allowed to stand for 30 min after the addition of water (20 mL). After this time, the mixture was extracted with acetonitrile (50 mL) using a homogenizer. The homogenate was filtered under reduced pressure through Celite filter aid, and the residue was re-homogenized with acetonitrile (20 mL) and filtered once more under reduced pressure. The volume of the combined extract was made up to 100 mL with acetonitrile. An aliquot of the extract (20 mL) was added to a 50 mL PTFE centrifuge tube containing sodium chloride (10 g) and phosphate buffer (20 mL, 0.5 mol L<sup>-1</sup>, pH 7.0). The mixture was shaken vigorously for 10 min and centrifuged at 3,000 rpm for 5 min. The resulting acetonitrile layer was loaded onto a preconditioned (10 mL acetonitrile) ODS

cartridge and eluted with an additional portion of acetonitrile (5 mL). The combined eluate was concentrated to approximately 1 mL using a rotary evaporator (<40°C) and evaporated to dryness under a stream of nitrogen. The residue was subsequently redissolved in acetonitrile (3 mL) assisted by ultrasonication and an aliquot of toluene (1 mL) was added. The resulting solution was loaded onto a GCB/PSA cartridge and eluted with a mixture of acetonitrile/toluene (18 mL, 3:1). The combined eluate was concentrated to approximately 1 mL using a rotary evaporator (<40°C), evaporated to dryness under a stream of nitrogen and the residue was redissolved in an acetone/hexane mixture (2 mL, 3:7). The resulting solution was loaded onto a silica gel cartridge and eluted with a mixture of acetone/hexane (13 mL, 3:7). The combined eluate was concentrated to approximately 1 mL using a rotary evaporator (<40°C), evaporated to dryness under a stream of nitrogen, and finally, the residue was redissolved in an acetone/hexane mixture (1 mL, 1:1).

#### Recovery test

The recovery tests were performed for green tea, oolong tea, black tea and matcha (powdered green tea), at 0.01 mg kg<sup>-1</sup> spike level, with each test being replicated five times for all matrices. An aliquot (0.5 mL) of the working mixed standard solution of 162 pesticides was added to each sample and the mixture was allowed to stand for 30 min before extraction. The recoveries of the individual pesticides were calculated on the basis of the peak areas of each pesticide, and they were quantified using matrix-matched calibration (25%, 50%, 75%, 100%, 125% and 150%).

#### Preparation of matrix-matched standards

Matrix-matched standards were prepared as follows. Extracts of blank samples (100  $\mu$ L) were evaporated to dryness under a stream of nitrogen, and the residue was redissolved in the standard acetone/hexane solution (100  $\mu$ L, 1:1).

## Results and discussion

#### Optimization of sample preparation

**Extraction and salting-out.** As acetonitrile is the most commonly used extraction solvent for the multiresidue method for pesticides and it is also used in the Japanese official multiresidue method (Fig. 1), it was chosen as the extraction solvent for our studies. The extraction efficiency of the incurred pesticide residues obtained from tea samples was compared by Cajka et al., both with and without matrix swelling prior to extraction with acetonitrile.<sup>16]</sup> They demonstrated that the addition of water to the sample matrix is a key factor for achieving the maximum

extraction yield. Therefore, in this study, we mixed water (20 mL) with tea leaves (5.0 g) and allowed them to swell for 30 min prior to extraction with acetonitrile. Kanrar et al. assessed the extractability of both polar and nonpolar pesticide residues by shaking, blending and vortexing prepared tea and spent leaves and comparing the results obtained.<sup>17</sup> Their studies revealed that blending gave better recoveries for the majority of pesticides, compared to shaking or vortex-based methods. Thus, in this study, we adopted a homogenizing procedure using acetonitrile for extraction, as in the Japanese method. In addition, as high pesticides recoveries were obtained in the salting-out procedure conducted with the aid of sodium chloride and phosphate buffer, and both water and polar matrix components were efficiently removed from the extract, we also chose to adopt a salting-out step, as in the Japanese official multiresidue method.

**Cleanup.** In order to remove matrix components of low polarity, cleanup using an ODS cartridge was investigated. The acetonitrile layer obtained after salting-out was loaded directly onto an ODS cartridge and eluted with an additional 5 mL of acetonitrile. Although the yellow pigments from the tea extract eluted from the ODS cartridge, the majority of green pigments (chlorophylls) and other matrix components of low polarity remained on the cartridge. All pesticides examined in the study eluted from the

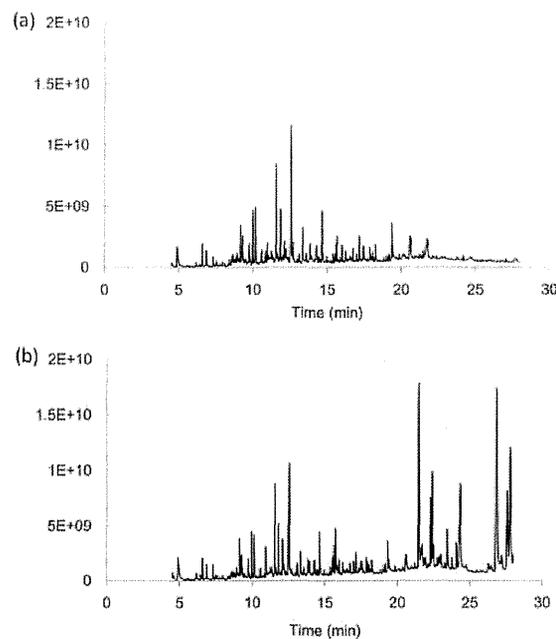


Fig. 2. Total ion current (TIC) chromatograms of the blank green tea sample, (a) with ODS cartridge cleanup and (b) without ODS cartridge cleanup.

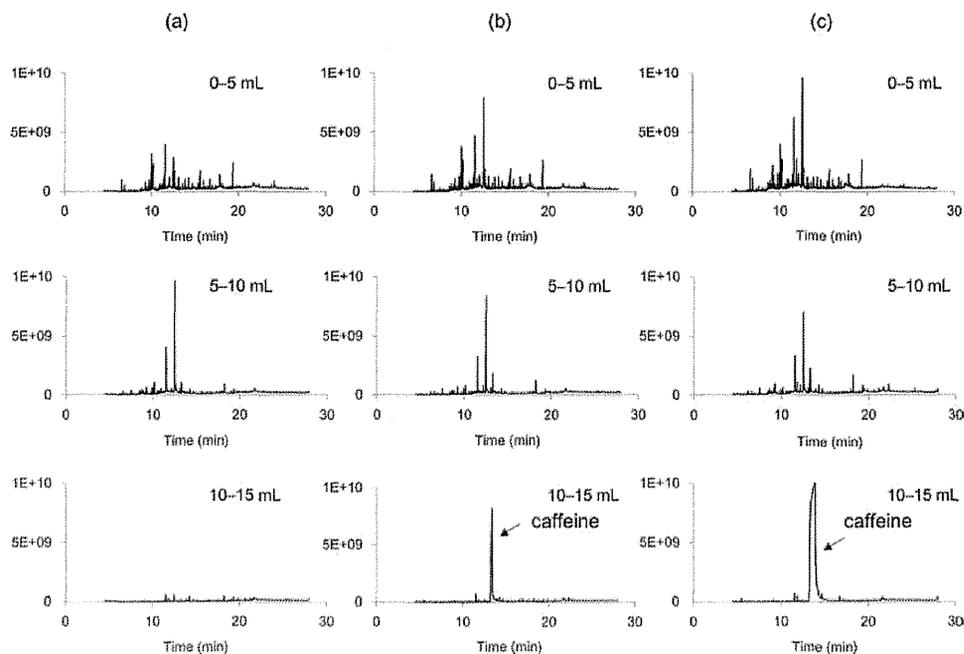


Fig. 3. Total ion current (TIC) chromatograms of a blank green tea sample after cleanup with a silica gel cartridge, eluting with (a) acetone/hexane (3:7), (b) acetone/hexane (4:6) and (c) acetone/hexane (5:5). Elution volume: upper 0–5 mL; middle 5–10 mL; lower 10–15 mL.

cartridge in high recoveries. Figure 2 shows the total ion current (TIC) chromatograms of green tea extract, both with and without ODS cartridge cleanup. It can be seen from Figure 2 that the number of coextracted matrix peaks (retention time >21 min) was smaller in the TIC chromatogram of the sample that had been subjected to ODS cartridge cleanup. Consequently, cleanup using an ODS cartridge was included in our method after the salting-out step, unlike in the Japanese official multiresidue method.

As shown in Figure 1, the Japanese official multiresidue method involves a tandem GCB/aminopropyl silanized silica gel (NH<sub>2</sub>) cartridge cleanup. However, as PSA is a stronger anion-exchange sorbent than the NH<sub>2</sub>-based sorbent, and can retain acidic coextracted matrix components

more efficiently, we selected the tandem GCB/PSA cartridge for use in this study.

In the cleanup procedure of the Japanese official multiresidue method (Fig. 1), the residue obtained from the salting-out step is dissolved in an aliquot of acetonitrile/toluene (2 mL, 3:1) to load onto the tandem GCB/NH<sub>2</sub> cartridge. However, because tea contains a large amount of matrix components, the residue obtained from the ODS cartridge cleanup was insoluble in the acetonitrile/toluene mixture. This was particularly the case for fermented black tea, most likely due to the presence of high-molecular-weight polyphenols, which could not be well partitioned in the aqueous phase during the salting-out step. In this study, we therefore chose to dissolve the residue in a small amount of a relatively high polar

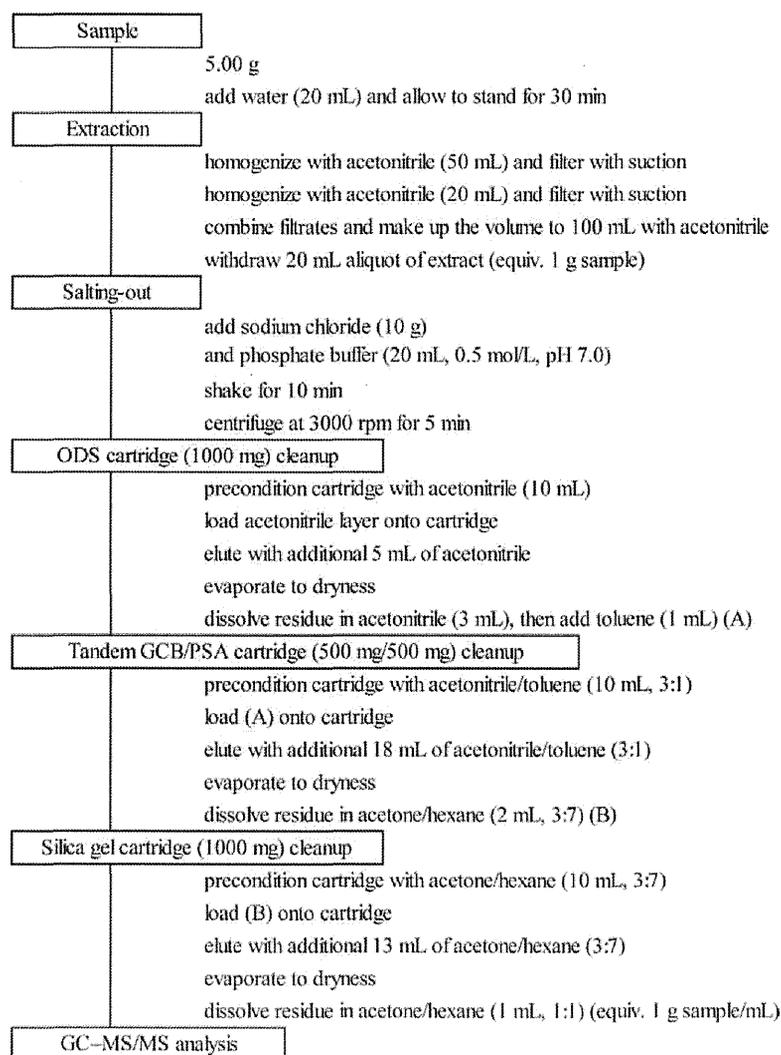


Fig. 4. Flow chart of sample preparation according to the modified multiresidue method.

solvent, acetonitrile (3 mL), with the aid of ultrasonication. The relatively apolar toluene (1 mL) was then added to the sample prior to loading onto the tandem GCB/PSA cartridge. While poor recoveries were obtained for coumaphos (50%), diflufenican (58%), pyrazophos (64%) and tecnazene (69%) using the standard mixture in pure solvent, all other pesticides were obtained in satisfactory recoveries from this system by elution with acetonitrile/toluene (3:1).

Although the tandem GCB/PSA cartridge cleanup effectively removed acidic matrix components together with other pigments from the samples, a large interfering peak with a retention time of 13.4 min was observed in the TIC chromatograms of tea extracts after the tandem GCB/PSA cartridge cleanup. This

peak is due to caffeine, which is present in large quantities (1–4%)<sup>[1]</sup> in tea extracts. The presence of such a significant quantity of caffeine in the final test solution for GC-MS/MS analyses would be troublesome, since this can affect the GC-MS/MS performance and cause considerable matrix effects. In order to avoid the requirement for frequent maintenance of GC-MS/MS systems, the introduction of a further cleanup step was investigated to reduce the caffeine content prior to GC-MS/MS analysis. As caffeine is a relatively polar compound compared to the pesticides being examined here, a silica gel cartridge was chosen for the separation of pesticides and caffeine, and the elution solvent was optimized by varying the proportion of acetone to hexane.

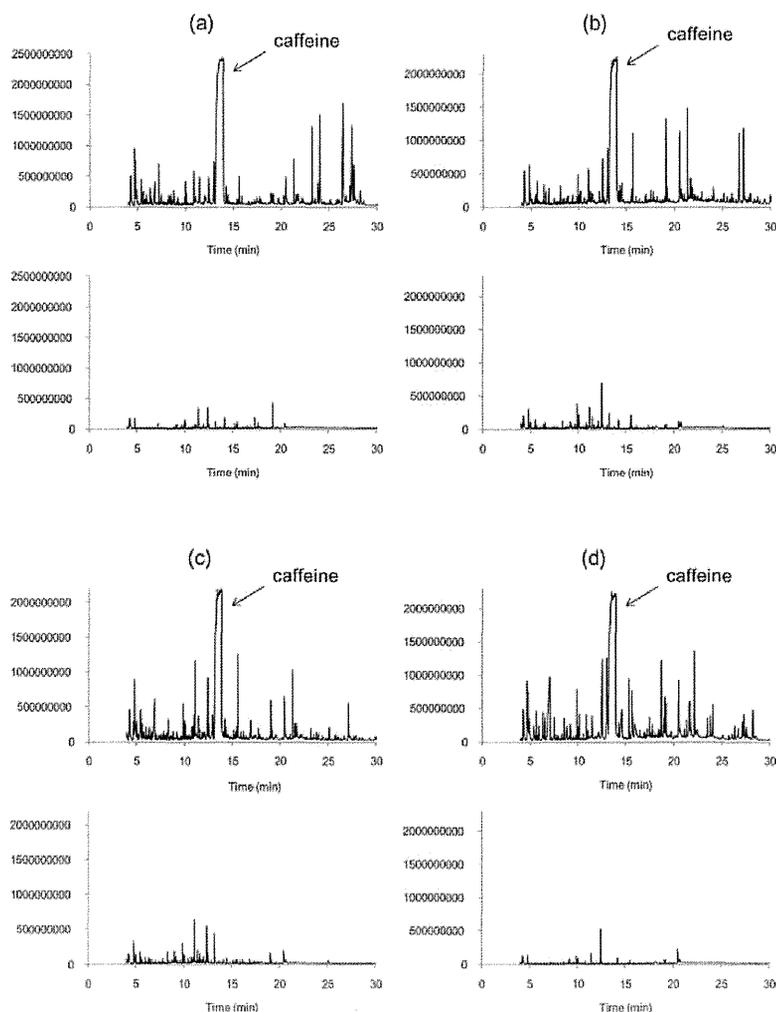


Fig. 5. Total ion current (TIC) chromatograms of (a) green tea, (b) oolong tea, (c) black tea and (d) matcha samples prepared by the Japanese official multiresidue method (upper) and the modified multiresidue method (lower).

Table 1. Recovery of pesticides from green tea, oolong tea, black tea and matcha (powdered green tea).

| Pesticide                  | Green tea    |      | Oolong tea   |      | Black tea    |      | Matcha (Powdered green tea) |      |
|----------------------------|--------------|------|--------------|------|--------------|------|-----------------------------|------|
|                            | Recovery (%) | RSD% | Recovery (%) | RSD% | Recovery (%) | RSD% | Recovery (%)                | RSD% |
| Acetochlor                 | 93           | 2    | 91           | 4    | 93           | 4    | 83                          | 4    |
| Acrinathrin                | 85           | 5    | 83           | 2    | 83           | 5    | 77                          | 3    |
| Alachlor                   | 99           | 2    | 92           | 3    | 91           | 3    | 83                          | 1    |
| Aldrin                     | 80           | 4    | 84           | 5    | 82           | 8    | 71                          | 2    |
| Ametryn                    | 96           | 7    | 89           | 5    | 87           | 3    | 80                          | 2    |
| Anilofos                   | 92           | 5    | 90           | 5    | 94           | 3    | 81                          | 7    |
| Aramite                    | 97           | 6    | 93           | 2    | 93           | 7    | 84                          | 3    |
| Atrazine                   | 100          | 7    | 92           | 2    | 90           | 8    | 82                          | 3    |
| Azinphos methyl            | 94           | 3    | 88           | 3    | 88           | 3    | 80                          | 3    |
| Azoxystrobin               | 93           | 4    | 90           | 5    | 90           | 3    | 82                          | 11   |
| Benalaxyl                  | 96           | 2    | 96           | 1    | 89           | 3    | 85                          | 4    |
| Benfluralin                | 86           | 6    | 89           | 2    | 93           | 3    | 82                          | 15   |
| Benfuresate                | 97           | 2    | 93           | 1    | 92           | 1    | 83                          | 6    |
| Benoxacor                  | 96           | 2    | 90           | 2    | 91           | 3    | 82                          | 4    |
| $\alpha$ -BHC              | 92           | 4    | 89           | 3    | 90           | 6    | 79                          | 4    |
| $\beta$ -BHC               | 100          | 5    | 88           | 2    | 88           | 4    | 81                          | 4    |
| $\gamma$ -BHC              | 95           | 5    | 91           | 3    | 89           | 5    | 80                          | 2    |
| $\delta$ -BHC              | 98           | 6    | 89           | 3    | 90           | 4    | 78                          | 3    |
| Bifenox                    | 93           | 5    | 91           | 5    | 82           | 4    | 83                          | 5    |
| Bifenthrin                 | 93           | 4    | 91           | 1    | 83           | 5    | 79                          | 2    |
| Bitertanol                 | 34           | 15   | 60           | 14   | 70           | 16   | 17                          | 4    |
| Bromobutide                | 95           | 1    | 91           | 3    | 90           | 5    | 95                          | 7    |
| Bromophos                  | 94           | 6    | 88           | 4    | 84           | 6    | 82                          | 13   |
| Bromophos ethyl            | 93           | 8    | 89           | 7    | 86           | 5    | 79                          | 16   |
| Bromopropylate             | 92           | 3    | 90           | 3    | 84           | 4    | 84                          | 11   |
| Bupirimate                 | 94           | 5    | 103          | 5    | 89           | 10   | 89                          | 7    |
| Butachlor                  | 92           | 2    | 94           | 3    | 87           | 5    | 80                          | 3    |
| Butafenacil                | 95           | 2    | 91           | 1    | 90           | 6    | 82                          | 9    |
| Butamifos                  | 90           | 2    | 92           | 5    | 92           | 2    | 81                          | 8    |
| Cadusafos                  | 95           | 1    | 91           | 2    | 88           | 5    | 88                          | 2    |
| Cafenstrole                | 88           | 3    | 88           | 2    | 88           | 5    | 79                          | 4    |
| Carfentrazone ethyl        | 96           | 3    | 92           | 2    | 95           | 3    | 93                          | 5    |
| Chlorbenside               | 92           | 4    | 84           | 2    | 85           | 1    | 73                          | 7    |
| Chlordane ( <i>cis</i> )   | 90           | 5    | 84           | 4    | 86           | 8    | 82                          | 13   |
| Chlordane ( <i>trans</i> ) | 89           | 6    | 81           | 2    | 82           | 9    | 82                          | 17   |
| Chlorfenson                | 97           | 6    | 92           | 3    | 92           | 5    | 84                          | 3    |
| Chlorfenvinphos            | 96           | 4    | 89           | 4    | 87           | 4    | 88                          | 12   |
| Chlorobenzilate            | 93           | 3    | 91           | 3    | 90           | 2    | 79                          | 6    |
| Chlorpropham               | 95           | 2    | 90           | 1    | 92           | 4    | 83                          | 4    |
| Chlorpyrifos               | 91           | 5    | 92           | 6    | 86           | 3    | 83                          | 13   |
| Chlorpyrifos methyl        | 92           | 3    | 89           | 5    | 92           | 4    | 81                          | 12   |
| Chlorthal dimethyl         | 93           | 10   | 87           | 5    | 87           | 3    | 84                          | 14   |
| Clomazone                  | 96           | 6    | 90           | 2    | 90           | 4    | 83                          | 3    |
| Clomeprop                  | 93           | 3    | 93           | 2    | 85           | 4    | 85                          | 9    |
| Coumaphos                  | 98           | 1    | 92           | 7    | 87           | 4    | 83                          | 4    |
| Cyanazine                  | 99           | 3    | 90           | 5    | 83           | 4    | 78                          | 6    |
| Cyflufenamid               | 94           | 16   | 87           | 9    | 91           | 12   | 78                          | 10   |
| Cyfluthrin                 | 93           | 3    | 90           | 2    | 91           | 4    | 83                          | 3    |
| Cyhalothrin                | 93           | 2    | 95           | 3    | 88           | 3    | 79                          | 5    |
| Cyproconazole              | 36           | 7    | 44           | 9    | 56           | 17   | 24                          | 3    |
| Cyprodinil                 | 92           | 1    | 89           | 2    | 79           | 8    | 79                          | 7    |
| Deltamethrin               | 86           | 4    | 87           | 3    | 88           | 3    | 79                          | 3    |
| Di-allate                  | 88           | 2    | 89           | 1    | 87           | 3    | 77                          | 5    |

(Continued on next page)

Table I. Recovery of pesticides from green tea, oolong tea, black tea and matcha (powdered green tea). (Continued)

| Pesticide            | Green tea    |      | Oolong tea   |      | Black tea    |      | Matcha (Powdered green tea) |      |
|----------------------|--------------|------|--------------|------|--------------|------|-----------------------------|------|
|                      | Recovery (%) | RSD% | Recovery (%) | RSD% | Recovery (%) | RSD% | Recovery (%)                | RSD% |
| Diazinon             | 95           | 5    | 101          | 3    | 92           | 8    | 88                          | 3    |
| Dichloran            | 90           | 5    | 89           | 2    | 87           | 4    | 80                          | 4    |
| Dieldrin             | 99           | 11   | 84           | 18   | 91           | 11   | 73                          | 6    |
| Disulfoton           | 95           | 5    | 90           | 8    | 90           | 4    | 86                          | 9    |
| Dimethametryn        | 96           | 3    | 93           | 1    | 88           | 5    | 83                          | 3    |
| Dimethenamid         | 96           | 3    | 93           | 1    | 88           | 3    | 82                          | 5    |
| Dimethoate           | 77           | 3    | 73           | 8    | 72           | 10   | 71                          | 2    |
| Dimethylvinphos (E)  | 95           | 4    | 90           | 3    | 90           | 4    | 83                          | 4    |
| Dimethylvinphos (Z)  | 93           | 2    | 91           | 2    | 89           | 5    | 80                          | 7    |
| Disulfoton           | 85           | 6    | 56           | 9    | 71           | 5    | 78                          | 5    |
| Dithiopyr            | 94           | 3    | 88           | 10   | 91           | 6    | 83                          | 10   |
| Edifenphos           | 90           | 4    | 90           | 7    | 87           | 3    | 80                          | 5    |
| $\alpha$ -Endosulfan | 94           | 8    | 85           | 10   | 88           | 9    | 80                          | 6    |
| $\beta$ -Endosulfan  | 100          | 6    | 90           | 4    | 93           | 12   | 81                          | 5    |
| Endosulfan sulfate   | 98           | 13   | 95           | 2    | 89           | 6    | 79                          | 4    |
| Endrin               | 92           | 6    | 86           | 10   | 93           | 7    | 74                          | 1    |
| EPN                  | 91           | 4    | 89           | 3    | 82           | 3    | 80                          | 3    |
| Esprocarb            | 95           | 2    | 91           | 2    | 92           | 6    | 83                          | 6    |
| Ethion               | 93           | 1    | 90           | 1    | 91           | 2    | 79                          | 6    |
| Ethoprophos          | 96           | 2    | 89           | 2    | 90           | 4    | 83                          | 2    |
| Etofenprox           | 93           | 5    | 82           | 3    | 88           | 5    | 79                          | 3    |
| Etoxazole            | 93           | 2    | 91           | 4    | 89           | 9    | 88                          | 12   |
| Fenarimol            | 112          | 7    | 92           | 2    | 89           | 4    | 80                          | 3    |
| Fenchlorphos         | 91           | 5    | 93           | 3    | 95           | 3    | 84                          | 13   |
| Fenitrothion         | 91           | 8    | 89           | 6    | 92           | 2    | 81                          | 9    |
| Fenoxanil            | 97           | 6    | 93           | 5    | 92           | 3    | 81                          | 9    |
| Fenpropimorph        | 76           | 7    | 69           | 9    | 66           | 12   | 59                          | 4    |
| Fenvalerate          | 90           | 4    | 90           | 1    | 88           | 4    | 80                          | 3    |
| Fipronil             | 92           | 2    | 94           | 7    | 87           | 5    | 81                          | 4    |
| Flamprop methyl      | 98           | 1    | 93           | 3    | 89           | 4    | 82                          | 7    |
| Flucythrinate        | 95           | 3    | 90           | 3    | 91           | 3    | 81                          | 4    |
| Fludioxonil          | 91           | 1    | 91           | 7    | 92           | 2    | 80                          | 10   |
| Fluquinconazole      | 95           | 1    | 91           | 4    | 85           | 4    | 81                          | 7    |
| Flutolanil           | 97           | 4    | 91           | 3    | 91           | 4    | 84                          | 7    |
| Fosthiazate          | 92           | 5    | 89           | 7    | 86           | 4    | 86                          | 5    |
| Fthalide             | 95           | 3    | 83           | 7    | 80           | 4    | 78                          | 8    |
| Indoxacarb           | 100          | 8    | 90           | 9    | 90           | 4    | 84                          | 3    |
| Iprobenfos           | 93           | 4    | 91           | 4    | 91           | 3    | 82                          | 1    |
| Isazophos            | 99           | 6    | 90           | 7    | 89           | 5    | 81                          | 5    |
| Isofenphos           | 96           | 2    | 92           | 2    | 92           | 4    | 83                          | 6    |
| Isofenphos oxon      | 91           | 6    | 92           | 6    | 88           | 5    | 84                          | 7    |
| Isoprocarb           | 86           | 3    | 89           | 2    | 89           | 1    | 81                          | 3    |
| Isoprothiolane       | 98           | 3    | 94           | 5    | 92           | 8    | 89                          | 9    |
| Isoxadifen ethyl     | 96           | 6    | 90           | 4    | 90           | 4    | 81                          | 8    |
| Isoxathion           | 108          | 3    | 92           | 8    | 84           | 6    | 88                          | 6    |
| Kresoxim methyl      | 99           | 1    | 90           | 3    | 91           | 4    | 82                          | 6    |
| Lenacil              | 89           | 5    | 84           | 5    | 83           | 4    | 76                          | 3    |
| Malathion            | 95           | 2    | 91           | 7    | 91           | 5    | 84                          | 2    |
| Mefenacet            | 98           | 4    | 90           | 2    | 89           | 5    | 81                          | 2    |
| Mefenpyr diethyl     | 96           | 2    | 92           | 4    | 91           | 1    | 82                          | 6    |
| Mepronil             | 92           | 4    | 91           | 5    | 92           | 2    | 85                          | 10   |
| Metalaxyl            | 91           | 1    | 87           | 7    | 84           | 4    | 78                          | 6    |
| Methidathion         | 91           | 3    | 91           | 2    | 91           | 5    | 83                          | 6    |

(Continued on next page)

Table 1. Recovery of pesticides from green tea, oolong tea, black tea and matcha (powdered green tea). (Continued)

| Pesticide              | Green tea    |      | Oolong tea   |      | Black tea    |      | Matcha (Powdered green tea) |      |
|------------------------|--------------|------|--------------|------|--------------|------|-----------------------------|------|
|                        | Recovery (%) | RSD% | Recovery (%) | RSD% | Recovery (%) | RSD% | Recovery (%)                | RSD% |
| Methoxychlor           | 106          | 4    | 91           | 9    | 85           | 4    | 84                          | 4    |
| Metolachlor            | 92           | 2    | 89           | 5    | 90           | 3    | 83                          | 2    |
| Myclobutanil           | 5            | 13   | 10           | 89   | 22           | 70   | 4                           | 0    |
| Oxadiazon              | 96           | 6    | 88           | 8    | 88           | 8    | 85                          | 6    |
| Oxadixyl               | 86           | 5    | 83           | 3    | 86           | 3    | 74                          | 7    |
| Pacllobutrazol         | 88           | 4    | 88           | 3    | 88           | 3    | 77                          | 6    |
| Parathion              | 90           | 2    | 92           | 12   | 86           | 8    | 83                          | 8    |
| Parathion methyl       | 94           | 3    | 92           | 9    | 91           | 5    | 84                          | 8    |
| Penconazole            | 32           | 24   | 67           | 7    | 72           | 12   | 13                          | 5    |
| Pendimethalin          | 91           | 4    | 89           | 4    | 86           | 4    | 76                          | 6    |
| Permethrin             | 93           | 6    | 88           | 6    | 88           | 5    | 81                          | 2    |
| Phenothrin             | 100          | 8    | 88           | 9    | 82           | 11   | 81                          | 3    |
| Phenthoate             | 88           | 3    | 90           | 2    | 93           | 3    | 80                          | 4    |
| Phosalone              | 94           | 3    | 88           | 4    | 89           | 2    | 82                          | 2    |
| Phosmet                | 91           | 5    | 83           | 5    | 77           | 5    | 76                          | 3    |
| Piperonyl butoxide     | 95           | 4    | 85           | 7    | 90           | 2    | 81                          | 4    |
| Procymidone            | 97           | 3    | 94           | 2    | 91           | 8    | 85                          | 11   |
| Profenofos             | 86           | 2    | 91           | 10   | 92           | 6    | 88                          | 11   |
| Prometryn              | 91           | 4    | 88           | 4    | 88           | 4    | 84                          | 3    |
| Propiconazole          | 91           | 3    | 90           | 8    | 86           | 6    | 68                          | 11   |
| Propoxur               | 85           | 3    | 86           | 5    | 86           | 4    | 81                          | 2    |
| Propyzamide            | 98           | 4    | 89           | 3    | 93           | 5    | 83                          | 3    |
| Prothiofos             | 89           | 6    | 87           | 6    | 86           | 2    | 82                          | 12   |
| Pyraflufen ethyl       | 93           | 8    | 93           | 12   | 88           | 3    | 79                          | 13   |
| Pyrazophos             | 94           | 3    | 90           | 6    | 90           | 2    | 80                          | 5    |
| Pyributicarb           | 92           | 4    | 87           | 9    | 89           | 7    | 81                          | 1    |
| Pyridaben              | 95           | 3    | 89           | 4    | 88           | 3    | 78                          | 11   |
| Pyridafenthion         | 91           | 3    | 90           | 2    | 88           | 7    | 78                          | 7    |
| Pyrifenox (E)          | 87           | 5    | 83           | 4    | 78           | 7    | 77                          | 5    |
| Pyrifenox (Z)          | 93           | 3    | 83           | 4    | 83           | 6    | 81                          | 5    |
| Pyrimethanil           | 95           | 7    | 86           | 4    | 84           | 4    | 79                          | 3    |
| Pyriminobac methyl (E) | 96           | 2    | 93           | 1    | 90           | 1    | 84                          | 5    |
| Pyriminobac methyl (Z) | 97           | 2    | 92           | 3    | 87           | 3    | 84                          | 7    |
| Pyriproxyfen           | 91           | 1    | 89           | 3    | 89           | 5    | 80                          | 3    |
| Quinoxifen             | 92           | 4    | 85           | 3    | 83           | 4    | 74                          | 8    |
| Quintozene             | 85           | 3    | 84           | 10   | 76           | 5    | 77                          | 8    |
| Simeconazole           | 12           | 14   | 31           | 41   | 51           | 33   | 1                           | 56   |
| Tebufenpyrad           | 97           | 2    | 92           | 8    | 93           | 4    | 89                          | 8    |
| Teenazene              | 80           | 2    | 82           | 4    | 81           | 1    | 76                          | 7    |
| Tefluthrin             | 95           | 5    | 86           | 4    | 89           | 5    | 81                          | 5    |
| Terbufos               | 90           | 2    | 84           | 3    | 85           | 5    | 79                          | 4    |
| Tetraclorvinphos       | 94           | 3    | 88           | 2    | 85           | 4    | 80                          | 4    |
| Tetradifon             | 97           | 6    | 89           | 6    | 86           | 6    | 81                          | 8    |
| Thenylchlor            | 93           | 3    | 92           | 3    | 89           | 4    | 87                          | 4    |
| Thiobencarb            | 97           | 3    | 86           | 7    | 93           | 6    | 77                          | 9    |
| Tolelofos methyl       | 94           | 4    | 90           | 2    | 90           | 3    | 78                          | 9    |
| Toifenpyrad            | 92           | 3    | 91           | 7    | 87           | 3    | 79                          | 3    |
| Triadimefon            | 89           | 9    | 95           | 5    | 90           | 5    | 83                          | 10   |
| Triadimenol            | 77           | 5    | 81           | 5    | 87           | 2    | 51                          | 18   |
| Tri-allate             | 91           | 2    | 89           | 4    | 87           | 2    | 77                          | 8    |
| Triazophos             | 92           | 6    | 94           | 4    | 92           | 3    | 80                          | 9    |
| Tribuphos              | 100          | 3    | 89           | 8    | 89           | 6    | 82                          | 6    |
| Trifloxystrobin        | 88           | 10   | 90           | 9    | 94           | 10   | 83                          | 6    |

(Continued on next page)

Table 1. Recovery of pesticides from green tea, oolong tea, black tea and matcha (powdered green tea). (Continued)

| Pesticide     | Green tea    |      | Oolong tea   |      | Black tea    |      | Matcha (Powdered green tea) |      |
|---------------|--------------|------|--------------|------|--------------|------|-----------------------------|------|
|               | Recovery (%) | RSD% | Recovery (%) | RSD% | Recovery (%) | RSD% | Recovery (%)                | RSD% |
| Trifluralin   | 87           | 6    | 93           | 7    | 88           | 6    | 80                          | 16   |
| Uniconazole P | 90           | 8    | 91           | 5    | 88           | 4    | 76                          | 8    |
| Vinclozolin   | 99           | 2    | 93           | 4    | 92           | 9    | 81                          | 7    |

Initially, pesticide recoveries obtained from the silica gel cartridge were compared using four different elution solvent mixtures, namely, acetone/hexane mixtures in the ratios of 2:8, 3:7, 4:6 and 5:5. The majority of pesticides were eluted using either 15 mL of the 3:7 mixture, 10 mL of the 4:6 mixture or 10 mL of the 5:5 mixture. In contrast, using an acetone/hexane ratio of 2:8, relatively low recoveries were obtained for a number of pesticides, even when 20 mL of the eluting solvent was utilized. The triazole-based pesticides (bitertanol, cyproconazole, myclobutanil, penconazole, simeconazole and triadimenol) showed low recoveries with all solvents.

With the optimized solvent mixtures in hand, we investigated the efficiency of the silica gel cartridge cleanup using green tea extract. The quantity of matrix components present in the green tea extract was evaluated by comparing the TIC chromatograms using different elution solvent mixtures. Figure 3 shows TIC chromatograms of the green tea extract after cleanup using a silica gel cartridge and eluting with (a) acetone/hexane (3:7), (b) acetone/hexane (4:6) and (c) acetone/hexane (5:5). As can be seen, caffeine did not elute with the acetone/hexane solution (15 mL, 3:7). However, using acetone/hexane (4:6) or acetone/hexane (5:5), caffeine was eluted in the 10–15 mL fractions. The TIC chromatograms of acetone/hexane (15 mL, 3:7), acetone/hexane (10 mL, 4:6) and acetone/hexane (10 mL, 5:5) were then compared. The results showed that acetone/hexane (15 mL, 3:7) appeared to contain the lowest quantities of matrix components, thus demonstrating that the silica gel cartridge cleanup eluting with 15 mL of acetone/hexane (3:7) yielded a cleaner final test solution than elution with 10 mL acetone/hexane (4:6) or 10 mL acetone/hexane (5:5).

From these results, acetone/hexane (15 mL, 3:7) was considered the optimal solvent mixture for silica gel cartridge cleanup, with the successful separation of pesticides and caffeine, elution of the majority of pesticides and retention of caffeine and other pigments on the cartridge. Out of all the pesticides tested, the triazole pesticides yielded the poorest recoveries, possibly due to low recoveries from the silica gel cartridge. However, quantification of the triazole pesticides is still possible, as they are LC-amenable and thus they can be quantified using a previously reported LC-MS/MS method<sup>[17]</sup> that did not require a silica gel cartridge cleanup step.

The overall schematic representation of the proposed modified method is shown in Figure 4. TIC chromatograms of green tea, oolong tea, black tea and matcha samples prepared by the Japanese official multiresidue method and the modified multiresidue method are shown in Figure 5 and the representative chromatograms of pesticides, where the retention times overlap with that of caffeine, can be seen in Figure A1. From the results presented in this study, we could conclude that the cleanup procedure for this method effectively removed coextracted matrix components such as pigments, caffeine and polyphenols, and gave a cleaner final test solution cleaner than that obtained by the original Japanese official multiresidue method. This final test solution was suitable for the subsequent GC-MS/MS analysis.

#### Recovery tests

The recovery tests were performed for green tea, oolong tea, black tea and matcha (powdered green tea) at a 0.01 mg kg<sup>-1</sup> spike level, with each test being replicated five times for all matrices. A 0.01 mg kg<sup>-1</sup> spike level was chosen for these tests, because a uniform limit of 0.01 mg kg<sup>-1</sup> is applied in Japan to all pesticides for which maximum residue limits (MRLs) are not established. The results of the recovery tests are summarized in Table 1. Out of 162 pesticides tested, the recoveries of 157 pesticides for green tea (97% of the tested pesticides), 155 for oolong tea (96%), 158 for black tea (98%) and 154 for matcha (95%) were within the acceptable range according to the criteria set by the Japanese guidelines.<sup>[18]</sup> These guidelines state that the recoveries of 70–120% and relative standard deviations less than 25% at the 0.01 mg kg<sup>-1</sup> spike level must be achieved. As mentioned previously, a possible reason for the low recoveries of triazole pesticides (i.e., bitertanol, cyproconazole, myclobutanil, penconazole, simeconazole and triadimenol) could be low recoveries from the silica gel cartridges. In the case of coumaphos, diflufenican, pyrazophos and teconazole, despite poor recoveries in pure solvents using the tandem GCB/PSA cartridge, high recoveries from tea itself were achieved, indicating that coexisting matrix components facilitated their elution from the cartridge. In addition, no interfering peaks were observed in the chromatograms of the blank