

samples at retention times similar to those of the pesticides, thus confirming the high selectivity of the developed method and its suitability for the quantitative analysis of the GC-amenable pesticide residues in tea.

Conclusions

In this study, by modifying the Japanese official multiresidue method, we developed an efficient and reliable multiresidue method for the determination of 162 pesticides in tea using GC-MS/MS, which is an analytical method that can determine pesticide contents without the need for continual and time-consuming maintenance. A combination of ODS, tandem GCB/PSA and silica gel cartridge cleanups effectively removed coextracted matrix components such as pigments, caffeine and polyphenols. The final test solutions obtained by the proposed method were suitable for the subsequent GC-MS/MS analysis. Recoveries and precisions of the majority of the 162 tested pesticides from the samples of green tea, oolong tea, black tea and matcha at a spiking level of 0.01 mg kg⁻¹ were within the acceptable range of the criteria of the Japanese guidelines. Therefore, the majority of the tested pesticides can be detected quantitatively by the modified method at 0.01 mg kg⁻¹, which is the uniform limit, established in Japan. Poor recoveries were obtained with triazole pesticides, possibly due to low recoveries from the silica gel cartridges. However, these triazole pesticides could be quantified using a previously reported LC-MS/MS method^[17] that does not utilize the silica gel cartridge cleanup procedure. Overall, the results obtained suggest that the developed method is suitable for the quantitative analysis of the GC-amenable pesticide residues in tea.

Funding

This work was supported by a Health and Labour Sciences Research Grant from the Ministry of Health, Labour and Welfare in Japan.

References

- [1] Chu, D.-C. Green tea—its cultivation, processing of the leaves for drinking materials, and kinds of green tea. In *Chemistry and Applications of Green Tea*, Yamamoto, T., Juneja, L. R., Chu, D.-C., Kim, M., Eds.; CRC Press: Boca Raton, FL, 1997; 1–11.
- [2] Yang, C.S.; Wang, X.; Lu, G.; Picinich, S.C. Cancer prevention by tea: animal studies, molecular mechanisms and human relevance. *Nat. Rev. Cancer* **2009**, *9*, 429–439.
- [3] Yang, C.S.; Chung, J.Y.; Yang, G.Y.; Chhabra, S.K.; Lee, M.J. Tea and tea polyphenols in cancer prevention. *J. Nutr.* **2000**, *130*, 472S–478S.
- [4] Sood, C.; Jaggi, S.; Kumar, V.; Ravindranath, S.D.; Shanker, A. How manufacturing processes affect the level of pesticide residues in tea. *J. Sci. Food Agric.* **2004**, *84*, 2123–2127.
- [5] Rajski, L.; Lozano, A.; Belmonte-Valles, N.; Uclés, A.; Uclés, S.; Mezcuza, M.; Fernández-Alba, A.R. Comparison of three multiresidue methods to analyse pesticides in green tea with liquid and gas chromatography/tandem mass spectrometry. *Analyst* **2013**, *138*, 921–931.
- [6] Cajka, T.; Sandy, C.; Bachanova, V.; Drabova, L.; Kalachova, K.; Pulkrabova, J.; Hajslova, J. Streamlining sample preparation and gas chromatography-tandem mass spectrometry analysis of multiple pesticide residues in tea. *Anal. Chim. Acta* **2012**, *743*, 51–60.
- [7] Kanrar, B.; Mandal, S.; Bhattacharyya, A. Validation and uncertainty analysis of a multiresidue method for 42 pesticides in made tea, tea infusion and spent leaves using ethyl acetate extraction and liquid chromatography-tandem mass spectrometry. *J. Chromatogr. A* **2010**, *1217*, 1926–1933.
- [8] Kanrar, B.; Mandal, S.; Bhattacharyya, A. Validation and uncertainty analysis of a multiresidue method for 67 pesticides in made tea, tea infusion, and spent leaves using ethyl acetate extraction and gas chromatography/mass spectrometry. *J. AOAC Int.* **2010**, *93*, 411–424.
- [9] Huang, Z.; Li, Y.; Chen, B.; Yao, S. Simultaneous determination of 102 pesticide residues in Chinese teas by gas chromatography-mass spectrometry. *J. Chromatogr. B* **2007**, *853*, 154–162.
- [10] Zhao, C.L.; Ding, R.Y.; Huo, L.G.; Li, H.D.; Dong, Z.; Wang, F. E.; Yang, G.S.; Lu, X.; Aboul-Enen, H.Y. Determination of pesticide residues in tea by gas chromatography/triple quadrupole mass spectrometry with solid-phase extraction. *J. AOAC Int.* **2014**, *97*, 1001–1006.
- [11] Pang, G.F.; Fan, C.L.; Chang, Q.Y.; Li, Y.; Kang, J.; Wang, W. W.; Cao, J.; Zhao, Y.B.; Li, N.; Li, Z.Y.; Chen, Z.M.; Luo, F.J.; Lou, Z.Y. High-throughput analytical techniques for multiresidue, multiclass determination of 653 pesticides and chemical pollutants in tea—Part III: Evaluation of the cleanup efficiency of an SPE cartridge newly developed for multiresidues in tea. *J. AOAC Int.* **2013**, *96*, 887–896.
- [12] Wang, J.; Chow, W.; Leung, D. Applications of LC/ESI-MS/MS and UHPLC/Qq-TOF-MS for the determination of 141 pesticides in tea. *J. AOAC Int.* **2011**, *94*, 1685–1714.
- [13] Steiniger, D.; Lu, G.; Butler, J.; Phillips, E.; Fintschenko, Y. Determination of multiresidue pesticides in green tea by using a modified QuEChERS extraction and ion-trap gas chromatography/mass spectrometry. *J. AOAC Int.* **2010**, *93*, 1169–1179.
- [14] Lozano, A.; Rajski, L.; Belmonte-Valles, N.; Uclés, A.; Uclés, S.; Mezcuza, M.; Fernández-Alba, A.R. Pesticide analysis in teas and chamomile by liquid chromatography and gas chromatography tandem mass spectrometry using a modified QuEChERS method: validation and pilot survey in real samples. *J. Chromatogr. A* **2012**, *1268*, 109–122.
- [15] Anastassiades, M.; Lehoutay, S.J.; Stajnbaher, D.; Schenck, F.J. Fast and easy multiresidue method employing acetonitrile extraction/partitioning and “dispersive solid-phase extraction” for the determination of pesticide residues in produce. *J. AOAC Int.* **2003**, *86*, 412–431.
- [16] Director Notice; Analytical methods for residual compositional substances of agricultural chemicals, feed additives, and veterinary drugs in food; Ministry of Health, Labour and Welfare, Tokyo, Japan; 2005; Syoku-An No. 1129002.
- [17] Saito, S.; Nemoto, S.; Matsuda, R. Simultaneous determination of pesticide residues in tea by LC-MS/MS—Modification of Japanese official multiresidue method. *Jpn. J. Food Chem. Safety* **2014**, *21*, 27–36.
- [18] Director Notice; Guideline for the validation of analytical methods for agricultural chemical residues in food; Ministry of Health, Labour and Welfare, Tokyo, Japan; 2010; Syoku-An No. 1224-1.

Appendix

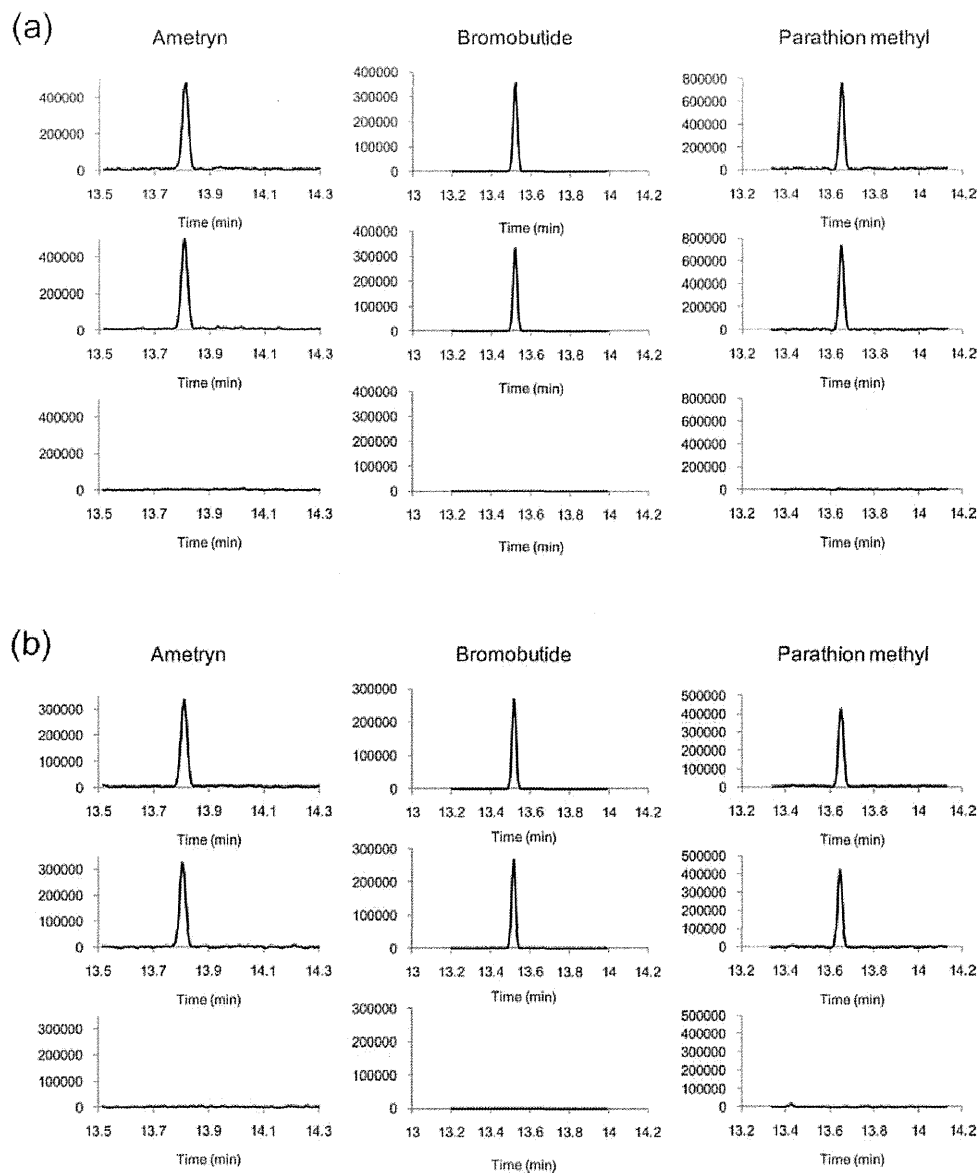


Fig. A1. Representative chromatograms of pesticides in (a) green tea, (b) oolong tea, (c) black tea and (d) matcha. Upper: matrix-matched standard. Middle: spiked sample. Lower: blank sample. (Continued)

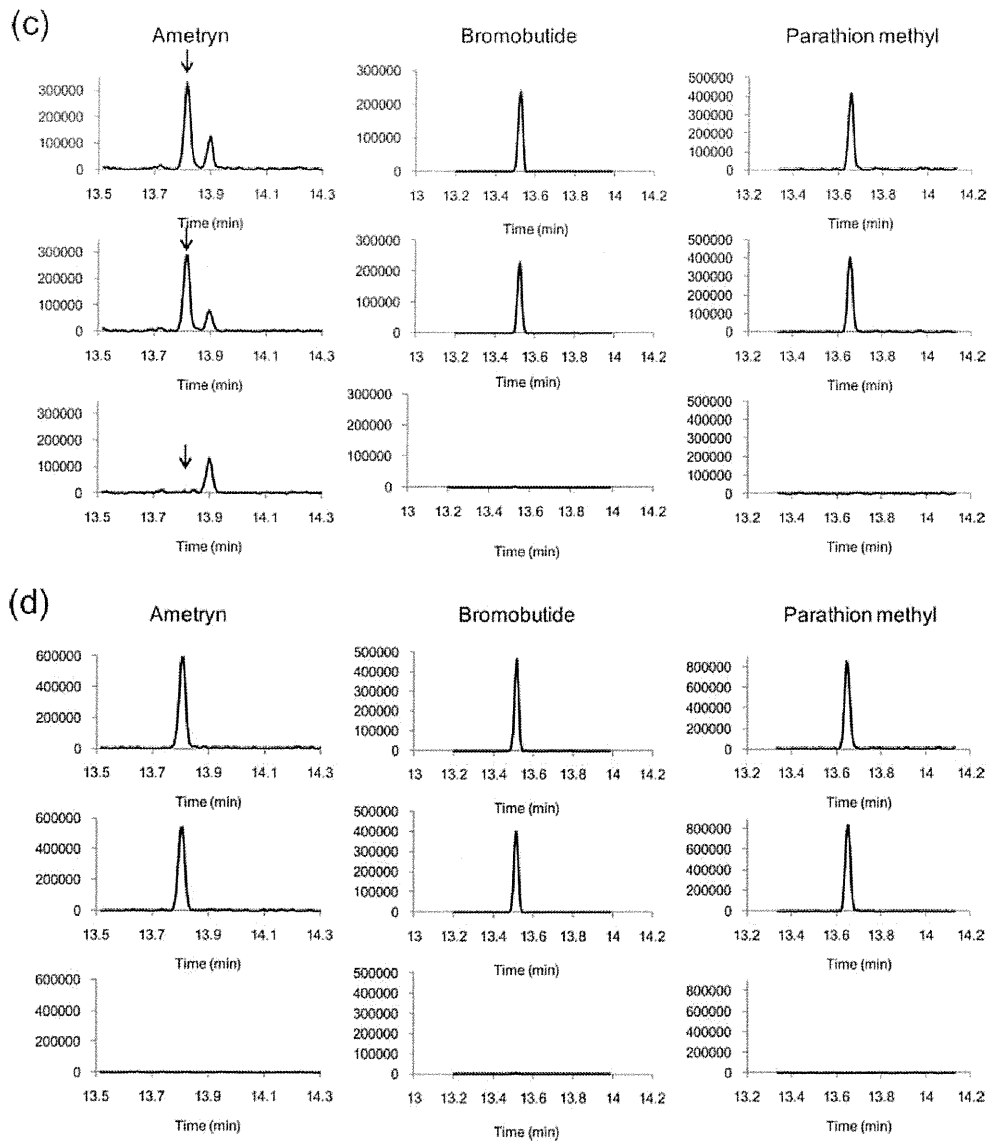


Fig. A1. Representative chromatograms of pesticides in (a) green tea, (b) oolong tea, (c) black tea and (d) matcha. Upper: matrix-matched standard. Middle: spiked sample. Lower: blank sample. (Continued from previous page)

Table A1. GC-MS/MS parameters for the tested pesticides.

Pesticide	Retention time (min)	Quantification ion			Confirmation ion		
		Precursor ion (m/z)	Product ion (m/z)	Collision energy (eV)	Precursor ion (m/z)	Product ion (m/z)	Collision energy (eV)
Acetochlor	13.51	223.1	132.1	14	223.1	147.1	14
Acerinathrin	19.94	208.0	181.1	8	181.0	152.1	20
Alachlor	13.67	188.1	160.1	8	237.1	160.1	15
Aldrin	14.60	262.8	193.0	24	262.8	228.0	20
Ametryn	13.83	227.1	170.1	10	227.1	185.1	10
Anilofos	19.17	226.1	157.0	20	226.1	184.1	8
Aramite	16.48, 16.54, 16.63, 16.83	185.0	63.0	14	334.1	185.1	8
Atrazine	12.15	200.1	122.1	10	215.1	200.2	10
Azinphos methyl	19.63	160.1	77.0	15	160.1	132.0	5
Azoxystrobin	23.71	344.1	329.3	18	388.1	360.3	8
Benalaxyl	17.71	234.1	174.1	10	266.1	148.1	10
Benfluralin	11.22	292.1	264.2	8	292.1	206.2	14
Benfuresate	13.39	256.1	163.1	8	163.1	121.0	5
Benoxacor	13.19	259.0	120.0	14	261.0	120.0	14
α -BHC	11.70	181.0	145.0	14	218.9	183.0	8
β -BHC	12.19	181.0	145.0	14	218.9	183.0	8
γ -BHC	12.96	218.9	183.0	8	181.0	145.0	14
δ -BHC	12.39	181.0	145.0	14	218.9	183.0	8
Bifenoxy	18.92	341.0	185.0	20	341.0	183.0	20
Bifenthrin	18.82	181.1	165.1	20	181.1	166.2	10
Bitertanol	20.60, 20.71	170.1	141.1	20	170.1	115.1	20
Bromobutide	13.53	232.2	114.1	8	232.2	176.2	8
Bromophos	14.92	330.9	316.1	14	328.9	313.9	20
Bromophos ethyl	15.71	358.9	303.1	20	302.8	285.0	14
Bromopropylate	18.92	340.9	183.0	20	185.0	157.0	14
Bupirimate	16.54	273.2	193.2	14	316.2	208.2	8
Butachlor	15.88	237.1	160.1	10	176.1	147.1	10
Butafenacil	20.92	331.0	180.0	18	331.0	124.0	14
Butamifos	16.05	286.1	202.1	14	286.1	258.2	8
Cadusafos	11.48	159.0	96.9	14	159.0	130.9	8
Cafenstrole	21.11	100.1	72.0	8	188.2	119.1	20
Carfentrazone ethyl	17.64	312.1	151.0	18	330.0	269.1	18
Chlorbenside	15.81	268.0	125.0	14	125.0	89.0	20
Chlordane (<i>cis</i>)	16.07	372.8	266.1	20	374.8	266.0	20
Chlordane (<i>trans</i>)	15.82	372.8	266.1	20	374.8	266.0	20
Chlorfenson	16.27	175.0	111.0	14	302.0	175.0	8
Chlorfenvinphos	15.10, 15.31	323.0	267.1	14	267.0	159.0	14
Chlorobenzilate	17.08	251.0	139.0	14	139.0	111.0	14
Chlorpropham	11.15	213.1	171.0	8	213.1	127.0	14
Chlorpyrifos	14.46	314.0	258.0	14	196.9	169.0	14
Chlorpyrifos methyl	13.53	286.0	271.1	14	288.0	273.1	14
Chlorthal dimethyl	14.56	331.9	301.1	20	300.9	273.1	14
Clomazone	12.25	204.1	107.0	20	125.0	89.0	14
Clomeprop	19.28	288.1	169.1	18	290.1	171.0	18
Coumaphos	20.79	226.0	163.1	18	362.1	226.1	14
Cyanazine	14.54	225.1	189.2	14	212.1	123.0	20
Cyflufenamid	16.73	412.1	295.2	12	294.1	203.1	20
Cyfluthrin	21.18, 21.28, 21.34, 21.38	163.0	127.0	8	206.1	151.0	20
Cyhalothrin	19.63, 19.80	197.1	141.0	14	181.1	152.1	14
Cyproconazole	16.88	222.1	125.0	20	224.1	127.0	20

(Continued on next page)

Table A1. GC-MS/MS parameters for the tested pesticides. (Continued)

Pesticide	Retention time (min)	Quantification ion			Confirmation ion		
		Precursor ion (<i>m/z</i>)	Product ion (<i>m/z</i>)	Collision energy (eV)	Precursor ion (<i>m/z</i>)	Product ion (<i>m/z</i>)	Collision energy (eV)
Cyprodimil	15.15	224.2	208.2	20	225.2	210.2	14
Deltamethrin	23.48	252.9	93.0	20	181.1	152.0	20
Di-allate	11.55, 11.75	234.1	150.0	20	234.1	192.1	14
Diazinon	12.54	179.1	137.1	20	199.1	135.1	14
Dichloran	11.94	206.0	176.0	14	176.0	148.0	20
Dieldrin	16.61	276.9	241.1	8	262.9	193.0	30
Diflufenican	18.23	394.1	266.2	18	266.1	246.2	14
Dimethametryn	15.24	212.1	122.1	12	212.1	94.1	20
Dimethenamid	13.42	230.1	154.1	8	203.1	126.0	20
Dimethoate	11.94	229.0	87.0	8	125.0	79.0	8
Dimethylvinphos (E)	14.23	295.0	109.0	16	297.0	109.0	16
Dimethylvinphos (Z)	14.51	295.0	109.0	16	297.0	109.0	16
Disulfoton	12.80	186.0	97.0	14	186.0	153.0	8
Dithiopyr	13.89	354.1	306.3	8	354.1	286.2	14
Edifenphos	17.84	310.1	109.0	20	173.0	109.0	14
α -Endosulfan	16.09	240.9	206.0	14	195.0	159.0	8
β -Endosulfan	17.22	240.9	206.0	14	195.0	159.0	8
Endosulfan sulfate	17.94	273.8	239.0	20	271.8	237.0	20
Endrin	17.02	262.9	228.0	18	262.9	193.0	30
EPN	18.89	157.0	110.0	14	169.1	141.0	8
Esprocarb	14.33	222.2	91.0	14	162.2	91.0	8
Ethion	17.24	231.0	128.9	20	231.0	175.0	14
Ethoprophos	10.92	158.0	96.9	16	200.0	158.0	8
Etofenprox	21.87	163.1	107.0	20	163.1	135.1	8
Etozazole	18.99	300.1	270.3	20	330.1	300.3	20
Fenarimol	20.09	251.0	139.0	20	219.1	107.0	12
Fenchlorphos	13.87	285.0	270.1	20	287.0	272.1	20
Fenitrothion	14.16	277.1	260.2	8	260.0	125.0	14
Fenoxanil	16.89	293.1	155.0	20	189.0	125.0	14
Fenpropimorph	14.60	303.3	128.1	8	128.1	110.1	14
Fenvalerate	22.51, 22.75	167.1	125.0	14	225.1	119.0	14
Fipronil	15.15	367.0	213.0	24	369.0	215.0	24
Flamprop methyl	16.47	230.1	170.1	14	276.1	105.0	8
Flucythrinate	21.68, 21.88	199.1	107.0	20	199.1	157.1	8
Fludioxonil	16.25	248.1	127.0	20	248.1	154.0	16
Fluquinconazole	20.81	340.1	298.0	14	340.1	313.0	14
Flutolanil	16.20	323.1	173.1	20	173.0	145.1	14
Fosthiazate	14.96, 15.02	195.0	103.0	8	195.0	139.0	8
Fthalide	14.86	242.9	215.0	14	271.9	243.0	14
Indoxacarb	23.23	203.0	134.0	14	218.1	203.1	14
Iprobenfos	13.10	204.0	91.0	14	246.1	204.1	8
Isazophos	12.80	285.0	161.1	14	257.0	119.0	20
Isofenphos	15.28	213.1	121.0	20	255.1	213.1	8
Isofenphos oxon	14.64	229.1	201.1	8	201.0	121.0	14
Isoprocarb	9.81	136.1	121.1	10	121.1	103.0	10
Isoprothiolane	16.32	290.1	118.0	14	231.1	189.1	8
Isoxadifen ethyl	17.66	204.0	176.1	20	222.1	205.2	8
Isoxathion	16.81	177.1	130.1	14	313.1	177.1	8
Kresoxim methyl	16.57	206.2	116.0	8	206.2	131.1	14
Lenacil	17.90	153.1	136.1	14	153.1	110.0	14
Malathion	14.32	173.1	127.0	8	127.1	99.0	8
Mefenacet	19.78	192.0	136.0	14	192.0	109.0	20

(Continued on next page)

Table A1. GC-MS/MS parameters for the tested pesticides. (Continued)

Pesticide	Retention time (min)	Quantification ion			Confirmation ion		
		Precursor ion (m/z)	Product ion (m/z)	Collision energy (eV)	Precursor ion (m/z)	Product ion (m/z)	Collision energy (eV)
Mefenpyr diethyl	18.49	253.0	189.1	20	255.0	189.1	20
Mepromil	17.49	269.2	119.1	14	269.2	210.2	8
Metaxyl	13.82	234.2	146.1	14	206.2	132.1	14
Methidathion	15.74	145.1	85.0	8	145.1	58.0	14
Methoxychlor	19.00	227.1	212.2	14	227.1	169.1	20
Metolachlor	14.43	238.1	162.2	14	162.1	133.1	14
Myclobutanil	16.53	179.1	125.0	20	288.2	179.1	14
Oxadiazon	16.42	258.0	175.0	8	175.0	112.0	14
Oxadixyl	17.22	233.2	146.1	8	163.1	132.1	14
Paclobutrazol	15.91	236.1	125.0	14	236.1	167.1	20
Parathion	14.63	291.1	109.0	14	291.1	136.1	8
Parathion methyl	13.67	263.1	109.0	14	263.1	124.9	20
Penconazole	15.25	248.1	157.1	20	159.0	123.0	14
Pendimethalin	15.11	252.1	162.1	14	252.1	191.2	8
Permethrin	20.65, 20.78	163.0	127.0	8	183.1	153.1	14
Phenothrin	19.29, 19.41	183.1	153.1	14	183.1	165.2	8
Phenthoate	15.44	274.0	121.0	12	274.0	246.2	8
Phosalone	19.52	182.0	111.0	14	182.0	138.1	8
Phosmet	18.86	160.1	77.0	20	160.1	133.0	14
Piperonyl butoxide	18.34	176.1	131.1	14	176.1	103.0	20
Procymidone	15.51	283.1	96.0	14	283.1	255.2	14
Profenofos	16.38	337.0	267.0	20	339.0	269.0	14
Prometryn	13.88	241.1	184.2	14	226.1	184.2	8
Propiconazole	17.82, 17.94	259.1	173.0	14	261.1	175.0	14
Propoxur	10.61	152.1	110.0	8	110.1	82.0	20
Propyzamide	12.50	173.0	145.0	14	175.0	147.0	14
Prothiofos	16.29	309.0	239.1	14	267.0	239.0	8
Pyraflufen ethyl	17.87	412.0	349.3	14	349.0	307.2	14
Pyrazophos	20.01	221.1	193.2	12	232.1	204.2	14
Pyributicarb	18.56	165.1	108.0	12	181.1	108.0	14
Pyridaben	20.84	309.1	147.2	20	364.2	147.1	20
Pyridafenthion	18.68	340.1	199.0	12	340.1	203.1	20
Pyrifenox (E)	15.83	262.0	227.2	8	262.0	200.1	14
Pyrifenox (Z)	15.32	262.0	200.0	20	262.0	192.0	20
Pyrimethanil	12.67	198.1	183.2	16	199.1	158.1	25
Pyriminobac methyl (E)	17.84	302.1	256.2	14	302.1	230.2	14
Pyriminobac methyl (Z)	17.02	302.1	256.2	14	302.1	230.2	14
Pyriproxyfen	19.68	136.1	96.0	14	136.1	78.0	20
Quinoxifen	17.86	237.1	208.2	20	272.0	237.2	8
Quintozene	12.26	236.9	142.9	20	294.9	237.0	20
Simeconazole	13.68	195.1	75.0	14	211.1	121.0	14
Tebufenpyrad	19.13	318.2	131.1	20	333.2	171.1	14
Tecnazene	10.47	260.9	203.0	14	202.9	142.9	20
Tefluthrin	12.80	177.1	127.0	20	197.1	141.1	8
Terbufos	12.45	231.0	128.9	20	231.0	175.0	14
Tetrachlorvinphos	15.83	328.9	109.0	20	330.9	109.0	20
Tetradifon	19.42	355.9	229.0	8	228.9	201.0	14
Thenylchlor	18.15	288.1	141.0	14	288.1	174.2	8
Thiobencarb	14.52	257.1	100.1	14	100.1	72.0	8
Tolclofos methyl	13.69	265.0	250.1	16	265.0	220.1	20
Tolfenpyrad	24.12	383.2	171.1	20	383.2	181.1	14
Triadimefon	14.68	208.1	181.0	8	208.1	127.0	14

(Continued on next page)

Table A1. GC-MS/MS parameters for the tested pesticides. (Continued)

Pesticide	Retention time (min)	Quantification ion			Confirmation ion		
		Precursor ion (m/z)	Product ion (m/z)	Collision energy (eV)	Precursor ion (m/z)	Product ion (m/z)	Collision energy (eV)
Triadimenol	15.57, 15.67	168.2	70.0	10	168.2	112.0	5
Tri-allate	12.95	268.1	184.0	14	270.1	186.0	20
Triazophos	17.53	257.1	162.1	8	285.1	162.1	20
Tribuphos	16.50	202.0	146.9	8	258.1	146.9	14
Trifloxystrobin	17.72	222.1	190.2	8	190.1	130.0	8
Trifluralin	11.16	306.1	264.2	8	306.1	206.2	14
Uniconazole P	16.47	234.1	165.0	10	234.1	137.1	15
Vinclozolin	13.59	212.1	172.0	14	285.0	212.1	14

